


NSO (Nitrogen/Sulfur/Oxygen) Pesticides Mix # 3 on ZB-MR1**Column:** Zebron™ ZB-MultiResidue™-1, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G016-11**Oven Profile:** 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 5 min**Carrier Gas:** Constant Flow Helium, 0.7 mL/min**Injection:** On-Column :1 0.5 µL @ 43°C**Detection:** Mass Selective (MSD) (335°C)**Analyst Note:** Analytes 10 ppm in Acetone**Products used in this application:**

NSO (Nitrogen/Sulfur/Oxygen) Pesticides Mix # 3 on ZB-MR1

ANALYTES:

- 1 Vernolate
- 2 Pebulate
- 3 Tetrahydrophthalimide
- 4 Molinate
- 5 Phenmedipham
- 6 Promecarb
- 7 Terbumeton
- 8 Othilinone
- 9 Metribuzin
- 10 Nitrothal-isopropyl
- 11 Pendimethalin
- 12 Allethrin
- 13 Allethrin
- 14 Kresoxim-methyl
- 15 Tetramethrin
- 16 Tetramethrin
- 17 Fenazaquin

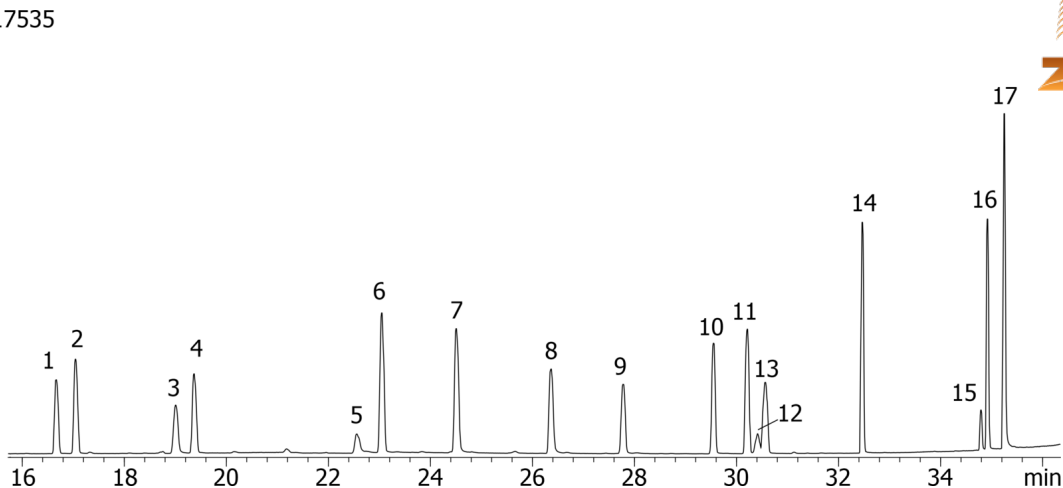