


**Organophosphorous Pesticides (EPA 8141) Mix # 2 on ZB-MR1****Column:** Zebron<sup>™</sup> ZB-MultiResidue<sup>™</sup>-1, GC Cap. Column 30 m x 0.25 mm x 0.25  $\mu$ m, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.25 mm x 0.25  $\mu$ m**Order No:** 7HG-G016-11**Oven Profile:** 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 5 min**Carrier Gas:** Constant Flow Helium, 0.7 mL/min**Injection:** On-Column :1 0.5  $\mu$ L @ 43°C**Detection:** Mass Selective (MSD) (335°C)**Analyst Note:** Analytes 10 ppm in Acetone**Products used in this application:**

17543

**ANALYTES:**

- 1 TEPP
- 2 Sulfotepp
- 3 Monocrotophos
- 4 Dimethoate
- 5 Malathion
- 6 Parathion
- 7 EPN

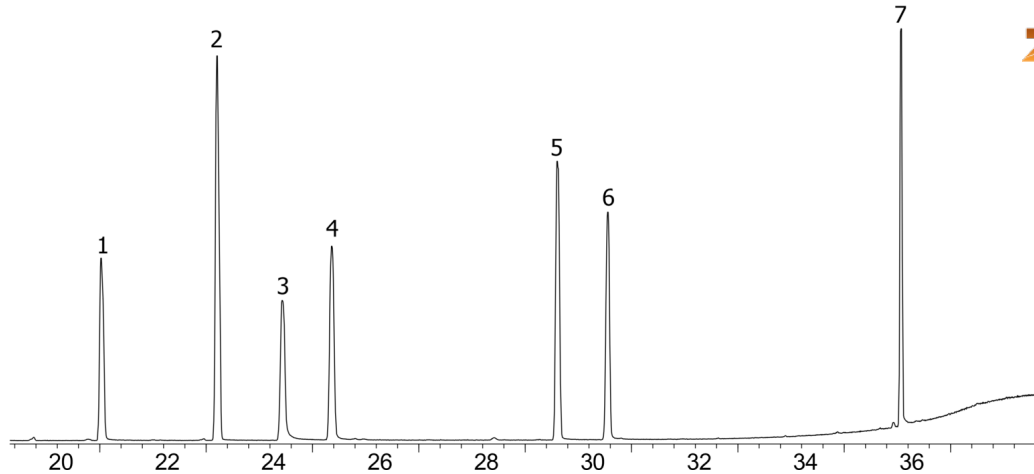