


European Organophosphorous Pesticides Mix # 1 on ZB-MR1**Column:** ZebronTM ZB-MultiResidueTM-1, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G016-11**Oven Profile:** 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 5 min**Carrier Gas:** Constant Flow Helium, 0.7 mL/min**Injection:** On-Column :1 0.5 μ L @ 43°C**Detection:** Mass Selective (MSD) (335°C)**Analyst Note:** Analytes 10 ppm in Acetone**Products used in this application:**

17546

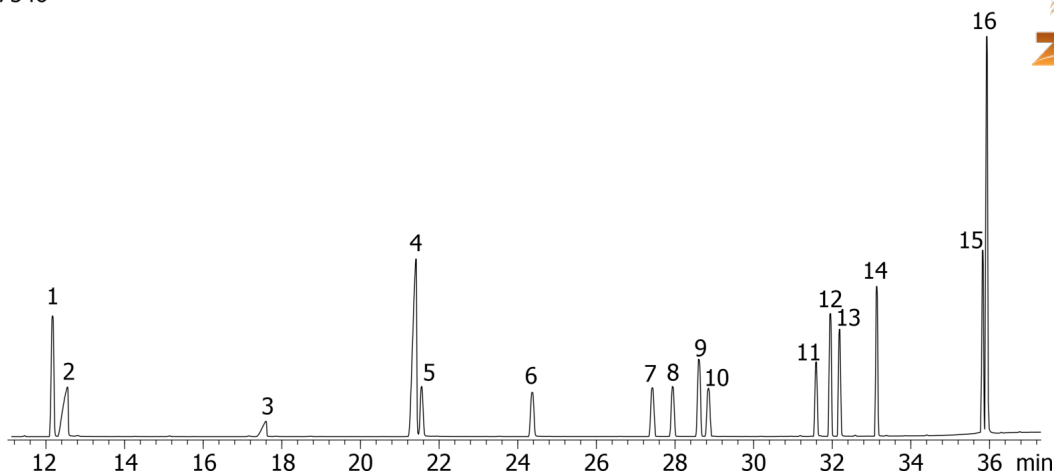

European Organophosphorous Pesticides Mix # 1 on ZB-MR1

ANALYTES:

- 1 Dichlorvos
- 2 Methamidophos
- 3 Acephate
- 4 Omethoate
- 5 Demeton-S-methyl
- 6 Dimethoate
- 7 Tolclofos-methyl
- 8 Pirimiphos methyl
- 9 Malathion
- 10 Chlorpyrifos
- 11 Methidathion
- 12 Prothiophos
- 13 Profenofos
- 14 Ethion
- 15 Azinphos-methyl
- 16 Pyrazophos

