


Organonitrogen Pesticides (EPA 633) Mix # 1 on ZB-MR1**Column:** Zebtron™ ZB-MultiResidue™-1, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G016-11**Oven Profile:** 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 5 min**Carrier Gas:** Constant Flow Helium, 0.7 mL/min**Injection:** On-Column :1 0.5 µL @ 43°C**Detection:** Mass Selective (MSD) (335°C)**Analyst Note:** Analytes 10 ppm in Acetone

17554

**Products used in this application:****ANALYTES:**

- 1 DEET
- 2 Terbacil
- 3 Metribuzin
- 4 Bromacil
- 5 Triadimefon
- 6 Tricyclazole
- 7 Hexazinone

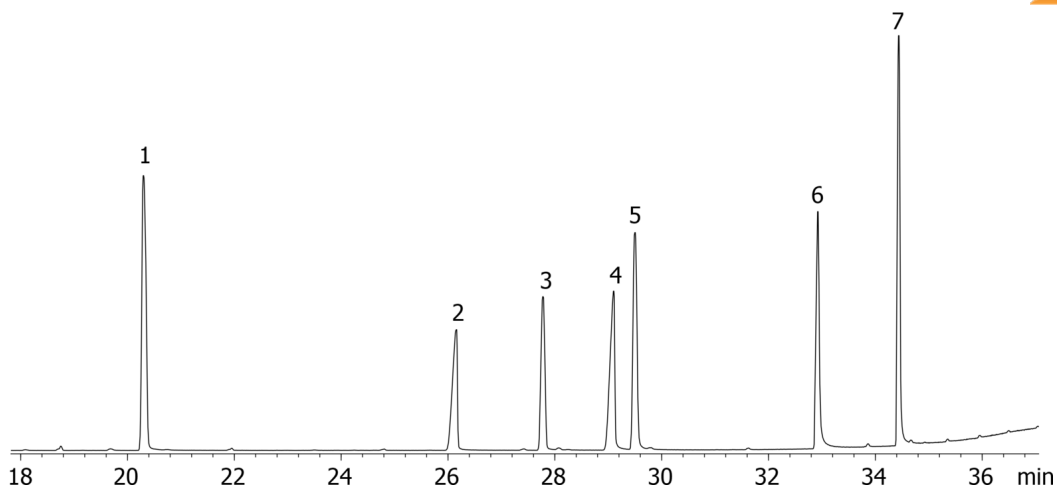