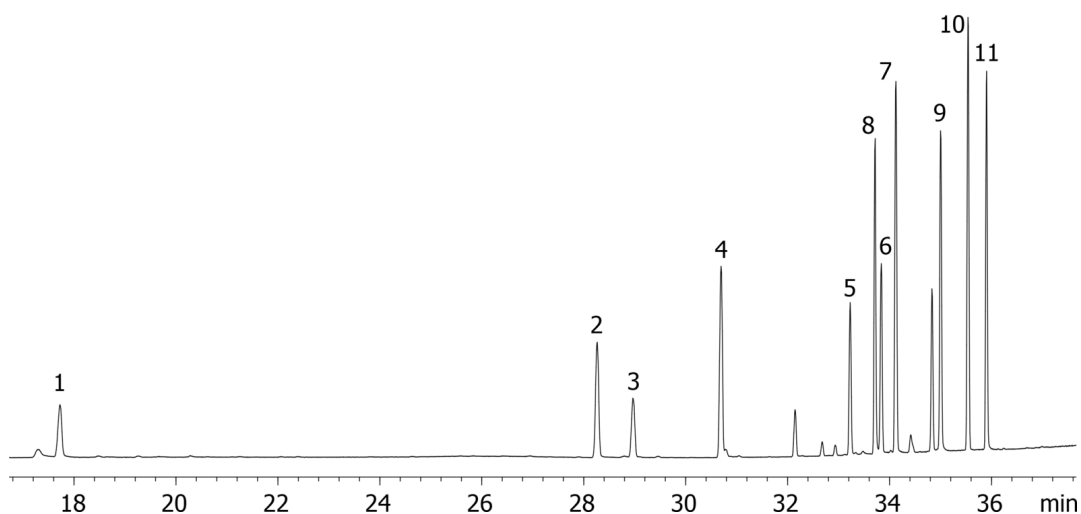


Organophosphorous Pesticides Mix 4 on ZB-MR2**Column:** Zebtron™ ZB-MultiResidue™-2, GC Cap. Column 30 m x 0.25 mm x 0.20 µm, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.25 mm x 0.2 µm**Order No:** 7HG-G017-10**Oven Profile:** 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 5 min**Carrier Gas:** Constant Flow Helium, 0.8 mL/min**Injection:** On-Column :1 0.5 µL @ 43°C**Detection:** Mass Selective (MSD) (335°C)**Analyst Note:** Analytes 10 ppm in Acetone

17565

**Products used in this application:**

Organophosphorous Pesticides Mix 4 on ZB-MR2

ANALYTES:

- 1 Trichlorfon
- 2 Parathion OA
- 3 Chlorpyrifos OA
- 4 Mecarbam
- 5 Chlorthiophos
- 6 Carbophenothion
- 7 Cyanofenphos
- 8 Tris(2-butoxyethyl) phosphate
- 9 Pyridaphenthion
- 10 Leptophos
- 11 Pyrazophos

