


Halogenated Pesticides Mix # 3 on ZB-MR2**Column:** Zebtron[™] ZB-MultiResidue[™]-2, GC Cap. Column 30 m x 0.25 mm x 0.20 μ m, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.25 mm x 0.2 μ m**Order No:** 7HG-G017-10**Oven Profile:** 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 5 min**Carrier Gas:** Constant Flow Helium, 0.8 mL/min**Injection:** On-Column :1 0.5 μ L @ 43°C**Detection:** Mass Selective (MSD) (335°C)**Analyst Note:** Analytes 10 ppm in Acetone**Products used in this application:**

Halogenated Pesticides Mix # 3 on ZB-MR2

ANALYTES:

- 1 Pentachlorobenzene
- 2 Tecnazene
- 3 Trifluralin
- 4 Chloropropham
- 5 BHC-alpha
- 6 Hexachlorobenzene
- 7 Dicloran
- 8 Quintozene
- 9 Pronamide
- 10 BHC-beta
- 11 Pentachloroaniline
- 12 BHC-delta
- 13 Vinclozolin
- 14 Dichlofluanid
- 15 Metolachlor
- 16 Chlorpyrifos
- 17 Dacthal
- 18 Tolyfluanid
- 19 Dieldrin
- 20 Trifloxystrobin
- 21 Quinoxifen
- 22 p,p'-DDT
- 23 Boscalid

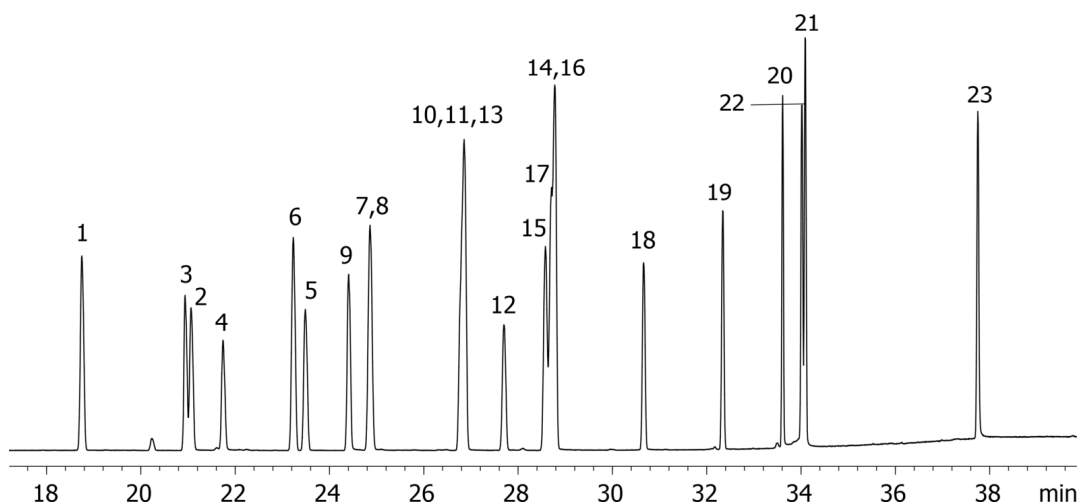