


NSO (Nitrogen/Sulfur/Oxygen) Pesticides Mix # 2 on ZB-MR2**Column:** Zebtron[™] ZB-MultiResidue[™]-2, GC Cap. Column 30 m x 0.25 mm x 0.20 μ m, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.25 mm x 0.2 μ m**Order No:** 7HG-G017-10**Oven Profile:** 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 5 min**Carrier Gas:** Constant Flow Helium, 0.8 mL/min**Injection:** On-Column :1 0.5 μ L @ 43°C**Detection:** Mass Selective (MSD) (335°C)**Analyst Note:** Analytes 10 ppm in Acetone**Products used in this application:**

NSO (Nitrogen/Sulfur/Oxygen) Pesticides Mix # 2 on ZB-MR2

ANALYTES:

- 1 Ethiolate
- 2 EPTC
- 3 Butylate
- 4 Propoxur
- 5 Cycloate
- 6 Bendiocarb
- 7 Carbofuran
- 8 Aminocarb
- 9 Fenfuram
- 10 Spiroxamine
- 11 Cymiazole
- 12 Ametryn
- 13 Ethofumesate
- 14 Fenpropimorph
- 15 Butralin
- 16 Aniten/Flurecol-butyl
- 17 Hexazinone
- 18 d-Phenothrin
- 19 Etofenprox
- 20 Famoxadone

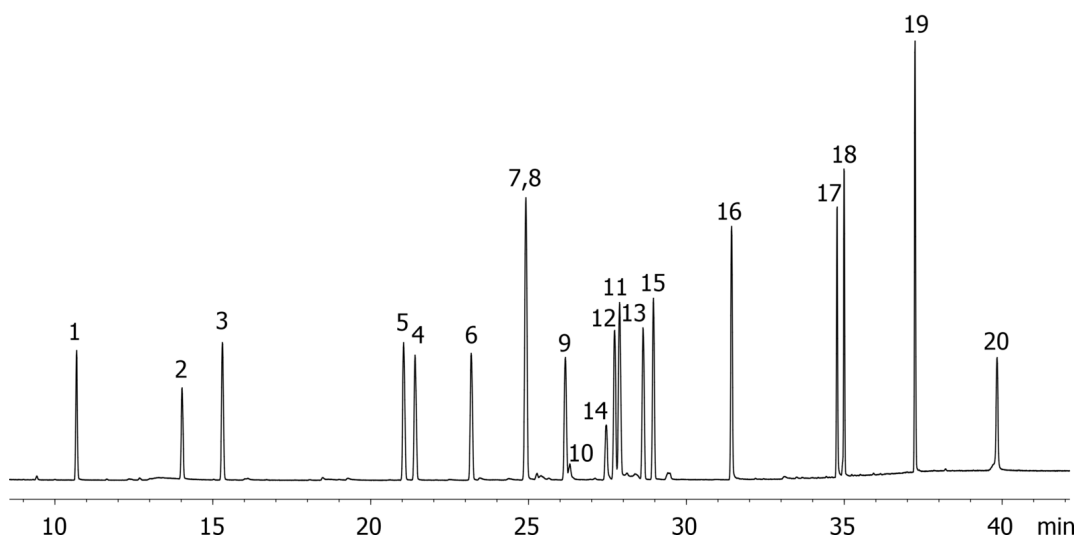