


Halogenated Pesticides Mix # 6 on ZB-MR2**Column:** Zebron[™] ZB-MultiResidue[™]-2, GC Cap. Column 30 m x 0.25 mm x 0.20 μ m, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.25 mm x 0.2 μ m**Order No:** 7HG-G017-10**Oven Profile:** 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 5 min**Carrier Gas:** Constant Flow Helium, 0.8 mL/min**Injection:** On-Column :1 0.5 μ L @ 43°C**Detection:** Mass Selective (MSD) (335°C)**Analyst Note:** Analytes 10 ppm in Acetone**Products used in this application:**

17573

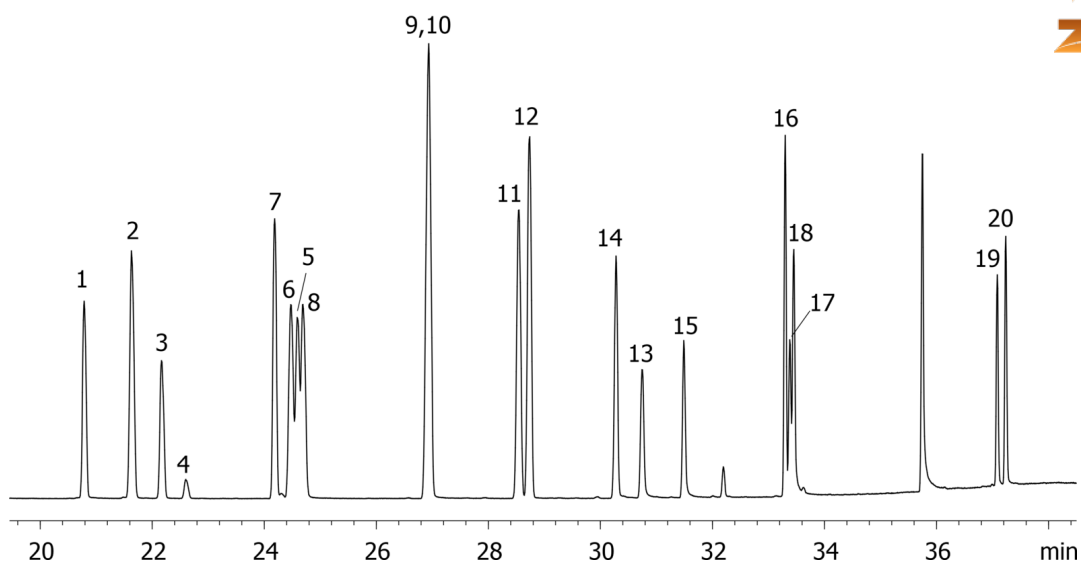

Halogenated Pesticides Mix # 6 on ZB-MR2

ANALYTES:

- 1 Ethalfluralin
- 2 2,3,5,6-Tetrachloroaniline
- 3 Di-allate
- 4 Di-allate
- 5 Simazine
- 6 Atrazine
- 7 Tefluthrin
- 8 Chlordene
- 9 Dimethachlor
- 10 Heptachlor
- 11 Pentachlorothioanisole
- 12 Dacthal
- 13 Pyrifenox
- 14 Fipronil
- 15 Pyrifenox
- 16 Nitrofen
- 17 Etaconazole
- 18 Etaconazole
- 19 Flucythrinate
- 20 Flucythrinate

