


Organophosphorous Pesticides (EPA 8141) Mix # 3 on ZB-MR2**Column:** Zebtron™ ZB-MultiResidue™-2, GC Cap. Column 30 m x 0.25 mm x 0.20 µm, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.25 mm x 0.2 µm**Order No:** 7HG-G017-10**Oven Profile:** 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 5 min**Carrier Gas:** Constant Flow Helium, 0.8 mL/min**Injection:** On-Column :1 0.5 µL @ 43°C**Detection:** Mass Selective (MSD) (335°C)**Analyst Note:** Analytes 10 ppm in Acetone**Products used in this application:**

17579

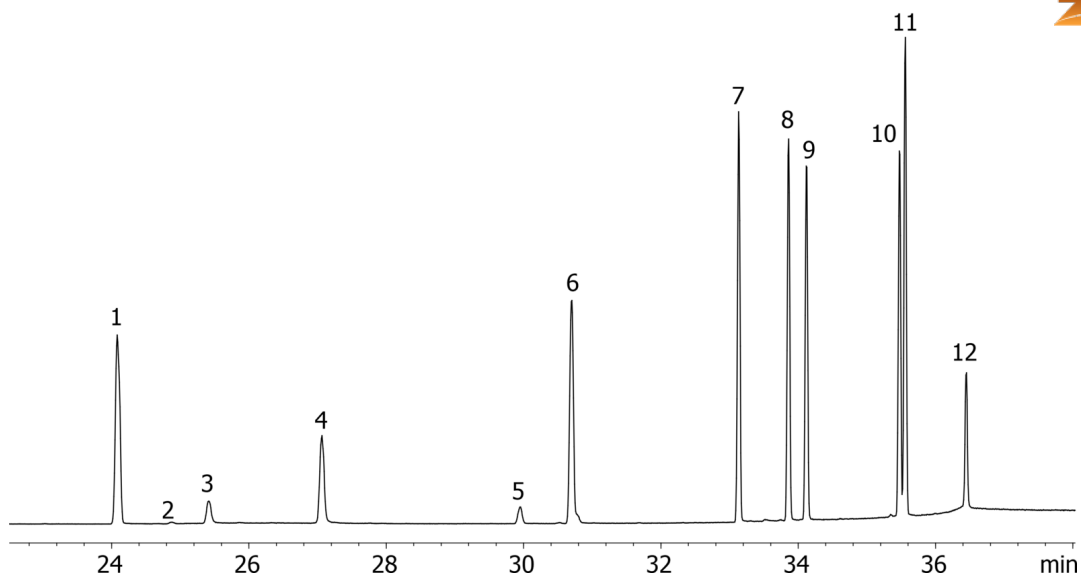

Organophosphorous Pesticides (EPA 8141) Mix # 3 on ZB-MR2

ANALYTES:

- 1 Terbufos
- 2 Dioxathion
- 3 Phosphamidon
- 4 Phosphamidon
- 5 Chlorfenvinphos
- 6 Chlorfenvinphos
- 7 Ethion
- 8 Carbophenothion
- 9 Famphur
- 10 Phosmet
- 11 Leptophos
- 12 Azinphos-ethyl

