


Organohalide Pesticides (EPA 505) Mix # 1 on ZB-MR2**Column:** Zebron™ ZB-MultiResidue™-2, GC Cap. Column 30 m x 0.25 mm x 0.20 µm, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.25 mm x 0.2 µm**Order No:** 7HG-G017-10**Oven Profile:** 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 5 min**Carrier Gas:** Constant Flow Helium, 0.8 mL/min**Injection:** On-Column :1 0.5 µL @ 43°C**Detection:** Mass Selective (MSD) (335°C)**Analyst Note:** Analytes 10 ppm in Acetone**Products used in this application:**

Organohalide Pesticides (EPA 505) Mix # 1 on ZB-MR2

ANALYTES:

- 1 Hexachlorocyclopentadiene
- 2 Hexachlorobenzene
- 3 Simazine
- 4 Atrazine
- 5 gama-BHC
- 6 Alachlor
- 7 Heptachlor
- 8 Aldrin
- 9 Heptachlor epoxide
- 10 Dieldrin
- 11 Endrin
- 12 Methoxychlor

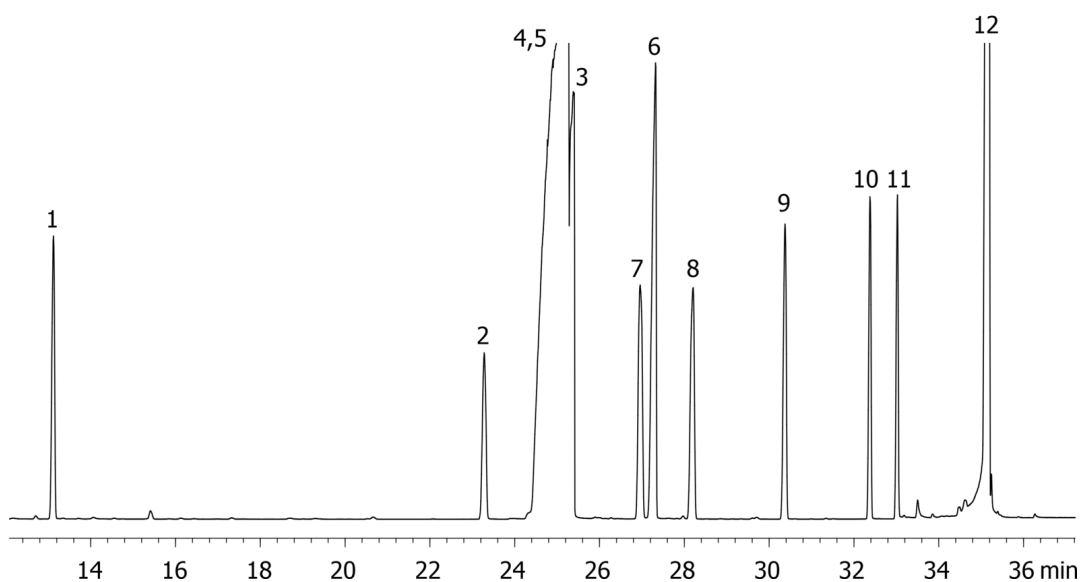