


European PAH Mix # 1 on ZB-MR2**Column:** Zebron[™] ZB-MultiResidue[™]-2, GC Cap. Column 30 m x 0.25 mm x 0.20 µm, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.25 mm x 0.2 µm**Order No:** 7HG-G017-10**Oven Profile:** 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 20 min**Carrier Gas:** Constant Flow Helium, 0.8 mL/min**Injection:** On-Column :1 0.5 µL @ 43°C**Detection:** Mass Selective (MSD) (335°C)**Analyst Note:** Analytes 10 ppm in Acetone**Products used in this application:**

17597

**ANALYTES:**

- 1 Benzo[j]fluoranthene
- 2 3-Methylcholanthrene
- 3 Dibenz[a,h]acridine
- 4 Dibenz[a,,j]acridine
- 5 7H-Dibenzo[c,g]carbazole
- 6 Dibenzo[a,e]pyrene
- 7 Dibenzo[a,i]pyrene
- 8 Dibenzo[a,h]pyrene

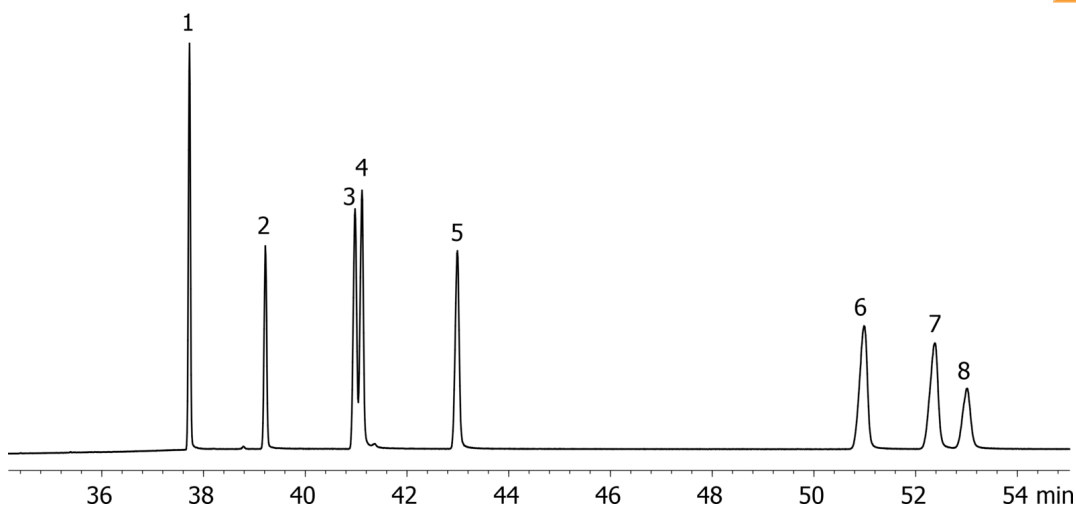