


Benzodiazepines by MSD on a ZB-XLB

Column: N/A
Phase: N/A
Dimensions: N/A
Order No: N/A
Oven Profile: N/A

Carrier Gas: N/A
Injection: N/A
Detection: N/A

Analyst Note: Analytes are 24ppm derivatized with MTBSTFA
Sim ions:
(Quant ion*) D5-Nordiazepam (*332,334): Nordiazepam (*327,328,383): D5-Oxazepam (*462,519): Oxazepam
(*457,513,459): D5-Temazepam (*362,288,390): Temazepam (*357,359,385): Lorazepam (*491,513,493):

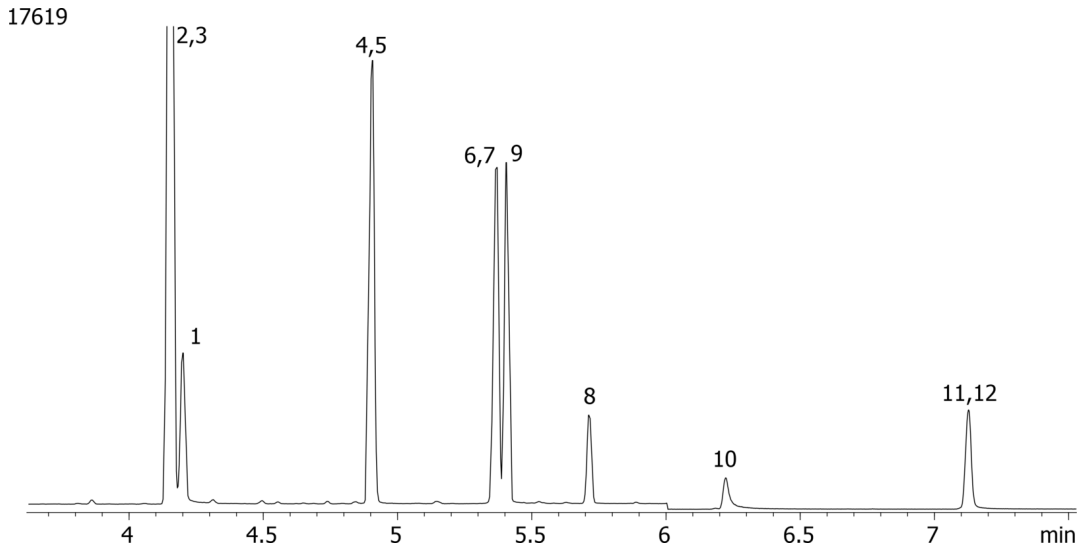

Benzodiazepines by MSD on a ZB-XLB

ANALYTES:

- 1 Diazepam
- 2 Nordiazepam-d5
- 3 Nordiazepam
- 4 Oxazepam-d5
- 5 Oxazepam
- 6 Temazepam-d5
- 7 Temazepam
- 8 Clonazepam
- 9 Lorazepam
- 10 Alprazolam
- 11 alpha-Hydroxyalprazolam-d5
- 12 alpha-Hydroxyalprazolam

