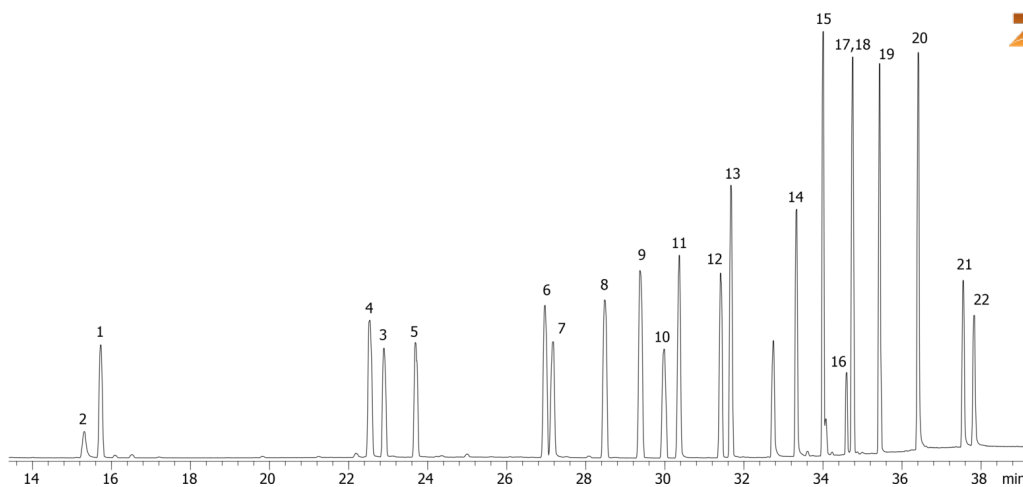


Halogenated Pesticides Mix # 5 on ZB-5ms

Column: Zebron[™] ZB-5MS, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea
Phase: 5% Phenyl-Arylene 95% Dimethylpolysiloxane
Dimensions: 30 meters x 0.25 mm x 0.25 µm
Order No: 7HG-G010-11
Oven Profile: 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 8 min
Carrier Gas: Constant Flow Helium, 0.8 mL/min
Injection: On-Column :1 0.5 µL @ 43°C
Detection: Mass Selective (MSD) (335°C)
Analyst Note: Analytes 10 ppm in Acetone

**Products used in this application:**

Halogenated Pesticides Mix # 5 on ZB-5ms

ANALYTES:

- 1 Etridiazole
- 2 3,4-Dichloroaniline
- 3 Profluralin
- 4 Clomazone
- 5 Dinitramine
- 6 Aldrin
- 7 Cyanazine
- 8 Heptachlor epoxide
- 9 Chlordane-trans
- 10 Nonachlor-trans
- 11 Flutolanil
- 12 Endrin
- 13 Chlorobenzilate
- 14 Epoxiconazole
- 15 Tebufenpyrad
- 16 Acrinathrin
- 17 Acrinathrin
- 18 Mirex
- 19 Fluquinconazole
- 20 Fluridone
- 21 Dimethomorph
- 22 Dimethomorph

