

HPLC Application

ID No.: 17708

Axia Packed Gemini 5 µm C18 (high pH loading - Lidocaine and Diphenhydramine - 22mg/1000µL)

Column: Gemini® 5 µm C18 110 Å, LC Column 50 x 21.2 mm, AXIA™ Packed, Ea

Dimensions: 50 x 21.2 mm ID

Order No: 00B-4435-P0-AX

Elution Type: Gradient

Eluent A: 0.2 % NH₄OH in Water

Eluent B: Acetonitrile

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	95	5
	2	5	5	95

Flow Rate: 30 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard PREP Cartridge, Gemini C18 15 x 21.2mm ID, Ea Part No.: AJ0-7846


- Holder Part No.: AJ0-8223



Products used in this application:



17708

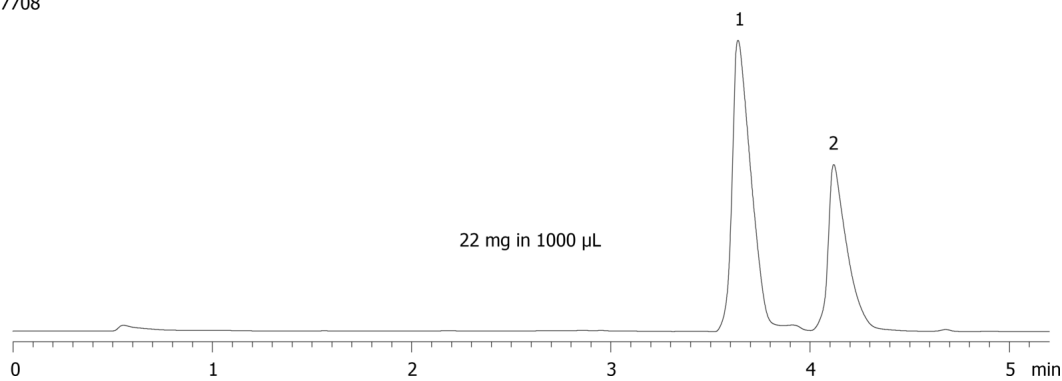

ANALYTES:

- 1 Lidocaine
- 2 Diphenhydramine

