

Loading Study: 60 mg at High pH on Gemini 5 µm C18 50x21.2 mm Axia Packed**Column:** Gemini® 5 µm C18 110 Å, LC Column 50 x 21.2 mm, AXIA™ Packed, Ea**Dimensions:** 50 x 21.2 mm ID**Order No:** 00B-4435-P0-AX**Elution Type:** Gradient**Eluent A:** 50 mM Ammonium Bicarbonate, pH 10.0**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	5	5	95

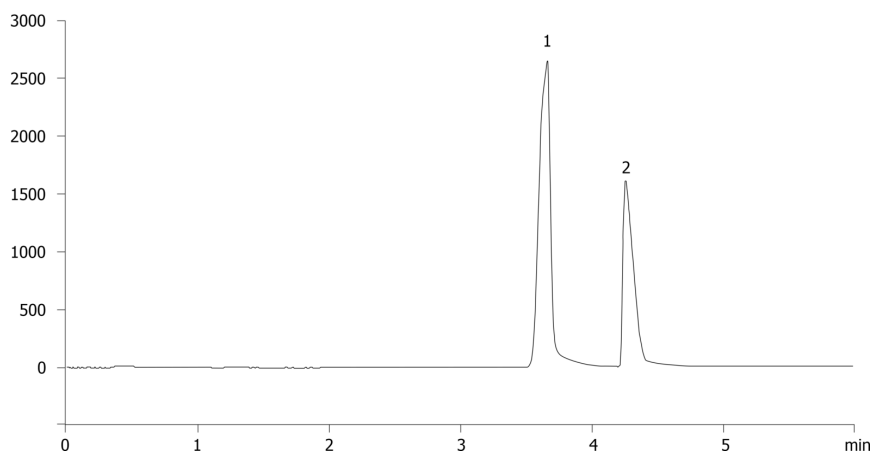
Flow Rate: 30 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard PREP Cartridge, Gemini C18 15 x 21.2mm ID, Ea Part No.: AJ0-7846

- Holder Part No.: AJ0-8223

**Products used in this application:**

17727

**ANALYTES:**

- 1 Propranolol
- 2 Diphenhydramine

