

Chlormezanone on Lux 5 μ Cellulose-1 250 x 4.6mm in PO-CH₃CN**Column:** Lux® 5 μ m Cellulose-1, LC Column 250 x 4.6 mm, Ea**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4459-E0**Elution Type:** Isocratic**Eluent A:** 0.1 % DEA in CH₃CN**Eluent B:** 0.1% DEA in IPA

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	10	95	5

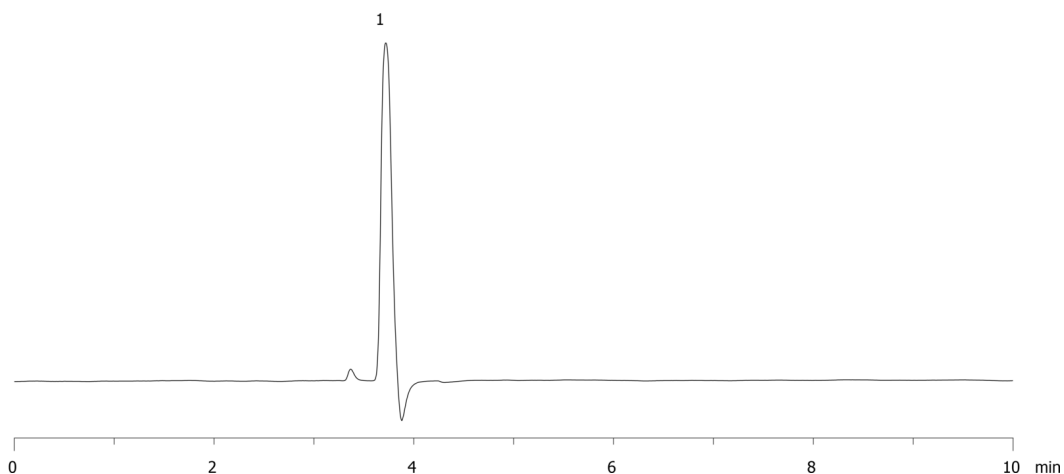
Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Cellulose-1, 4 x 3.0mm, 10/Pk Part No.: AJ0-8403

- Holder Part No.: KJ0-4282

**Products used in this application:**

17729

**ANALYTES:**

- 1** Chlormezanone
Retention Time: 3.7 min
- 2** Chlormezanone
Retention Time: 3.7 min

