


Propoxyphene by MSD on a ZB-5MSi**Column:** Zebtron™ ZB-5MSi, GC Cap. Column 10 m x 0.18 mm x 0.18 µm, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 10 meters x 0.18 mm x 0.18 µm**Order No:** 7CD-G018-08**Oven Profile:** 170 °C to 320 °C @ 20 °C/min**Carrier Gas:** Constant Flow Helium, 0.7 mL/min**Injection:** Split 10:1 1 µL @ 280°C**Detection:** Mass Selective (MSD) (320°C)**Analyst Note:** Analytes are 24ppm
Sim ions:
(Quant ion*) D11-Propoxyphene (*64,198,213): Propoxyphene (*58,193,208)**Products used in this application:**

17788

**ANALYTES:**

- 1 Propoxyphene-d11
- 2 Propoxyphene

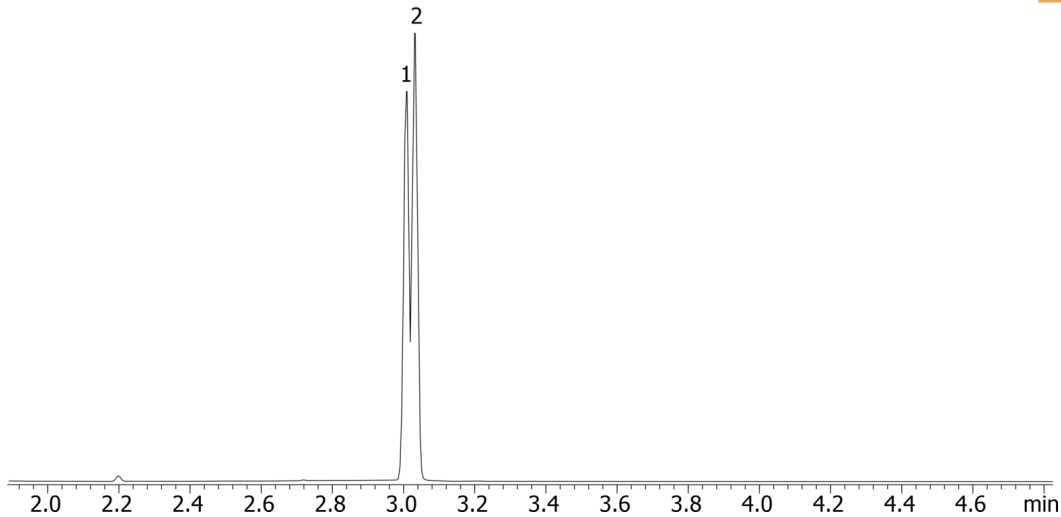