

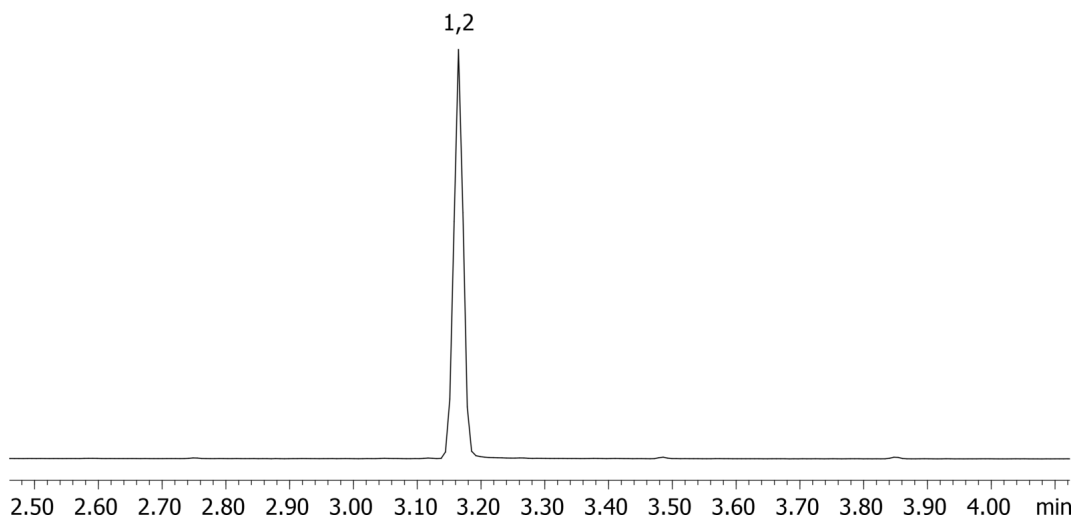
Benzoylecgonine by MSD on a ZB-Drug1

Column: Zebtron™ ZB-Drug-1, GC Cap. Column 10 m x 0.18 mm x 0.18 µm, Ea
Phase: Proprietary Drug Phase
Dimensions: 10 meters x 0.18 mm x 0.18 µm
Order No: 7CD-G023-08
Oven Profile: 235°C to 320°C @ 25 °C/min

Carrier Gas: Constant Flow Helium, 55 cm/sec
Injection: Split 10:1 1 µL @ 260°C
Detection: Mass Selective (MSD) (230°C)
Analyst Note: Analytes are 24ppm derivatized with BSTFA
Sim ions:
(Quant ion*) D3-Benzoylecgonine (*243,259,364): Benzoylecgonine (*240,256,361)

**Products used in this application:**

17805

**ANALYTES:**

- 1 Benzoylecgonine-d3
- 2 Benzoylecgonine

