


Opiates by MSD on a ZB-Drug1**Column:** Zebtron™ ZB-Drug-1, GC Cap. Column 10 m x 0.18 mm x 0.18 µm, Ea**Phase:****Dimensions:** 10 meters x 0.18 mm x 0.18 µm**Order No:** 7CD-G023-08**Oven Profile:** 180°C to 340°C @ 20 °C/min**Carrier Gas:** Constant Flow Helium, 55 cm/sec**Injection:** Split 10:1 1 µL @ 280°C**Detection:** Mass Selective (MSD) (230°C)**Analyst Note:** Analytes are 24ppm derivatized with BSTFA
Sim ions:
(Quant ion*) D6-Codeine (*377,349): Codeine (*371,343,234): D3-Morphine (*432,417): Morphine (*429,414,401)**Products used in this application:**

17806

**ANALYTES:**

- 1 Codeine-d6
- 2 Codeine
- 3 Morphine-d3
- 4 Morphine

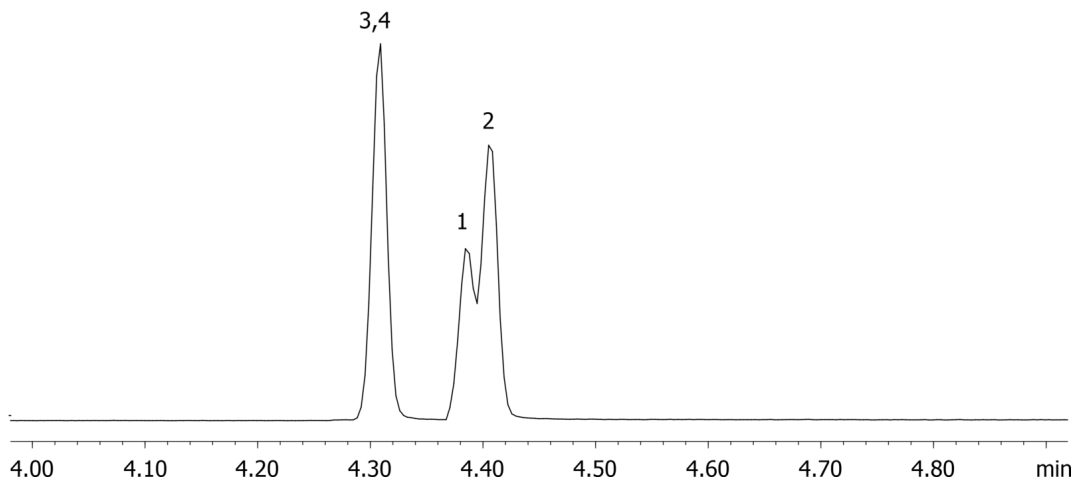