


Halogenated Pesticides Mix # 5 on ZB-50

Column: Zebron[™] ZB-50, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea
Phase: 50% Phenyl 50% Dimethylpolysiloxane
Dimensions: 30 meters x 0.25 mm x 0.25 μ m
Order No: 7HG-G004-11
Oven Profile: 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 8 min
Carrier Gas: Constant Flow Helium, 0.7 mL/min
Injection: On-Column :1 0.5 μ L @ 43°C
Detection: Mass Selective (MSD) (335°C)
Analyst Note: Analytes are 10 ppm in Acetone

**Products used in this application:**

17824

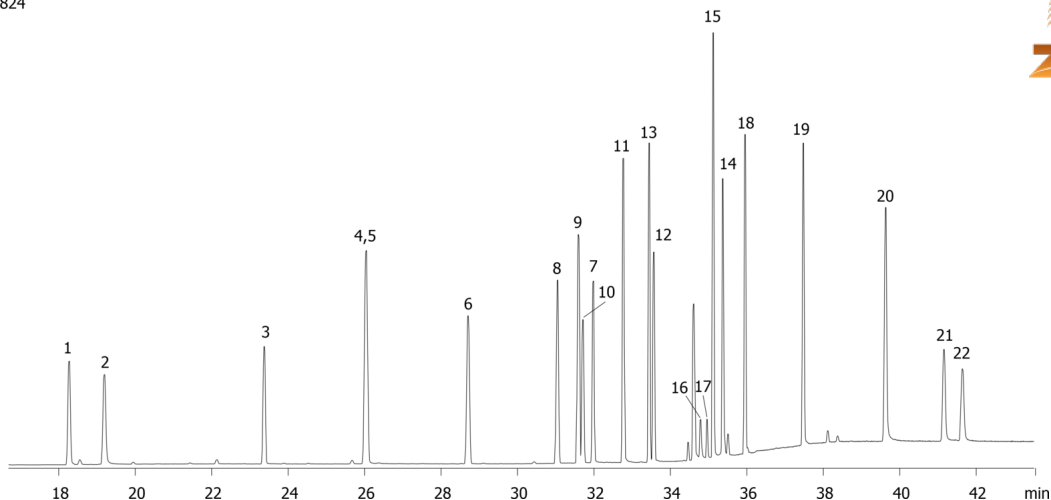

Halogenated Pesticides Mix # 5 on ZB-50

ANALYTES:

- 1 Etridiazole
- 2 3,4-Dichloroaniline
- 3 Profluralin
- 4 Clomazone
- 5 Dinitramine
- 6 Aldrin
- 7 Cyanazine
- 8 Heptachlor epoxide
- 9 Chlordane-trans
- 10 Nonachlor-trans
- 11 Flutolanil
- 12 Endrin
- 13 Chlorobenzilate
- 14 Epoxiconazole
- 15 Tebufenpyrad
- 16 Acrinathrin
- 17 Acrinathrin
- 18 Mirex
- 19 Fluquinconazole
- 20 Fluridone
- 21 Dimethomorph
- 22 Dimethomorph

