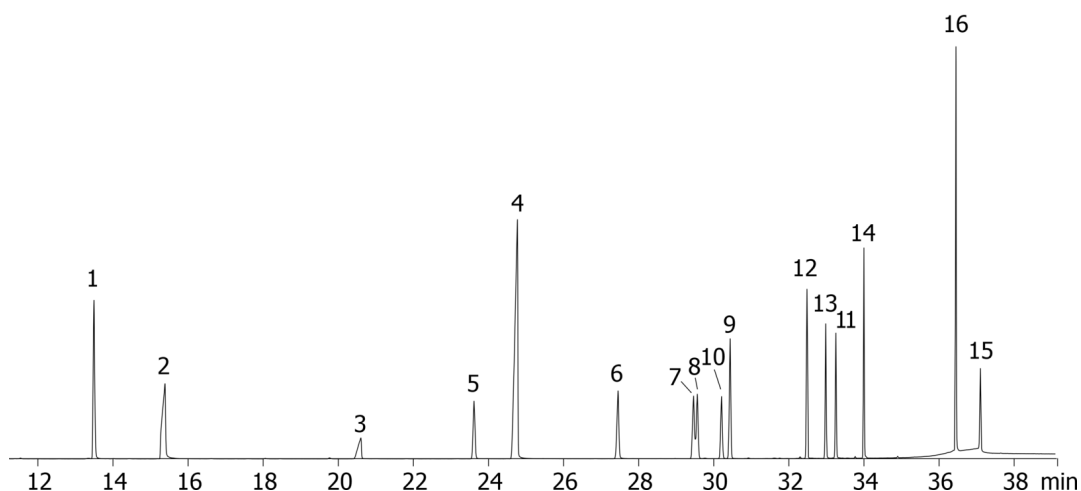


European Organophosphorous Pesticides Mix # 1 on ZB-50

Column: Zebron[™] ZB-50, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea
Phase: 50% Phenyl 50% Dimethylpolysiloxane
Dimensions: 30 meters x 0.25 mm x 0.25 μ m
Order No: 7HG-G004-11
Oven Profile: 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 8 min
Carrier Gas: Constant Flow Helium, 0.7 mL/min
Injection: On-Column :1 0.5 μ L @ 43°C
Detection: Mass Selective (MSD) (335°C)
Analyst Note: Analytes are 10 ppm in Acetone

17833

**Products used in this application:**

European Organophosphorous Pesticides Mix # 1 on ZB-50

ANALYTES:

- 1 Dichlorvos
- 2 Methamidophos
- 3 Acephate
- 4 Omethoate
- 5 Demeton-S-methyl
- 6 Dimethoate
- 7 Tolclofos-methyl
- 8 Pirimiphos methyl
- 9 Malathion
- 10 Chlorpyrifos
- 11 Methidathion
- 12 Prothiophos
- 13 Profenofos
- 14 Ethion
- 15 Azinphos-methyl
- 16 Pyrazophos

