

Pindolol on Lux 5 μ Cellulose-2 in RP 20% CH₃CN**Column:** Lux® 5 μ m Cellulose-2, LC Column 250 x 4.6 mm, Ea**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4457-E0**Elution Type:** Isocratic**Eluent A:** 20mM Ammonium acetate + 0.1%DEA**Eluent B:** Acetonitrile

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	10	80	20

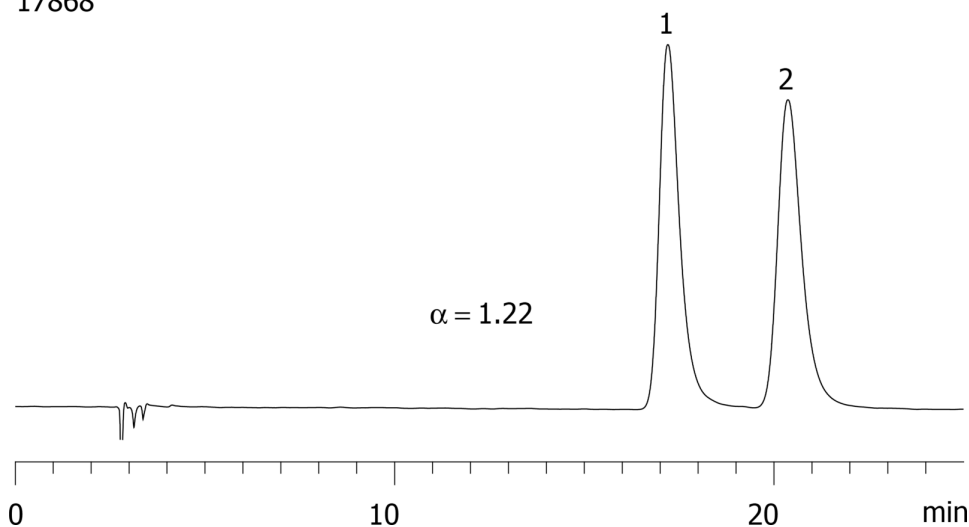
Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Cellulose-2, 4 x 3.0mm, 10/Pk Part No.: AJ0-8366

- Holder Part No.: KJ0-4282

**Products used in this application:**

17868

**ANALYTES:****1** Pindolol

Retention Time: 17.2 min

2 Pindolol

Retention Time: 20.4 min

