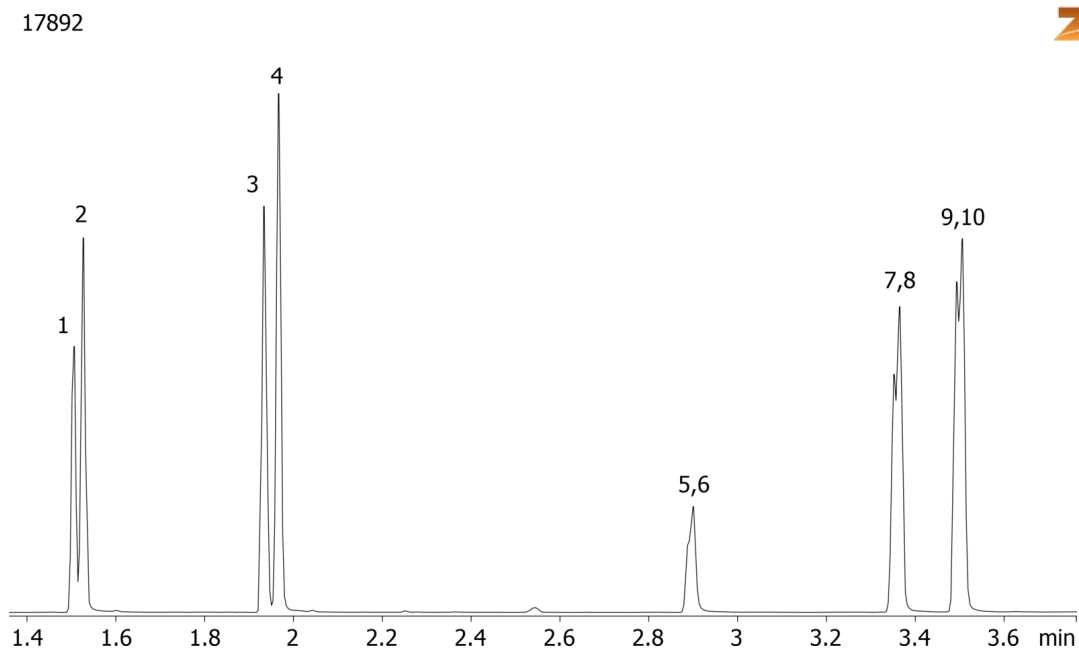


Expanded amphetamines by MSD on a ZB-5MS**Column:** Zebtron[™] ZB-5MS, GC Cap. Column 10 m x 0.18 mm x 0.18 μ m, Ea**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane**Dimensions:** 10 meters x 0.18 mm x 0.18 μ m**Order No:** 7CD-G010-08**Oven Profile:** 100 °C to 320 °C @ 25 °C/min**Carrier Gas:** Constant Flow Helium, 1 mL/min**Injection:** Split 50:1 0.5 μ L @ 280°C**Detection:** Mass Selective (MSD) (320°C)**Analyst Note:** Analytes are 40ppm derivatized with HFBA

Sim ions:

(Quant ion*) D11-Amphetamine (*244,128): Amphetamine (*240,118,117): D14-Methamphetamine (*261,213):
Methamphetamine (*254,210,118): D5-MDA (*167,380,244): MDA (*162,375,240): D5-MDMA (*258,213,164):**Products used in this application:**

Expanded amphetamines by MSD on a ZB-5MS

ANALYTES:

- 1 Amphetamine-d11
- 2 Amphetamine
- 3 Methamphetamine-d14
- 4 Methamphetamine
- 5 MDA-d5
- 6 MDA
- 7 MDMA-d5
- 8 MDMA
- 9 MDEA-d5
- 10 MDEA

