

Amphetamines in Urine by LC/MS/MS on Gemini 5µm C18 150 x 3.0

Column: Gemini® 5 µm C18 110 Å, LC Column 150 x 3 mm, Ea

Dimensions: 150 x 3 mm ID

Order No: 00F-4435-Y0

Elution Type: Isocratic

Eluent A: 10mM Ammonium formate w/0.1% formic acid

Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	0	90	10	

Flow Rate: 0.6 mL/min

Col. Temp.: ambient

Detection: Mass Spectrometer (MS) @ amu (450 °C)

Detector Info: <a target="_blank"

href="https://sciex.com/products/mass-spectrometers?utm_campaign=2019%20application%20search&utm_source=phenomenex&utm_medium=referral">SCIEX<

heater gas flow 8000cc/min
Heater temperature 450 °C; Positive Ion Mode


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Gemini C18 4 x 2.0mm, 10/Pk Part No.: AJ0-7596

- Holder Part No.: KJ0-4282



Products used in this application:

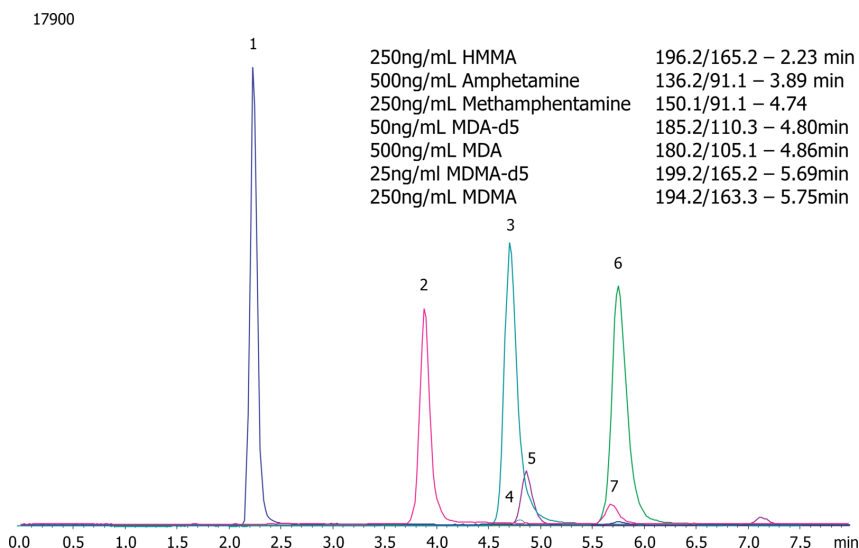

ANALYTES:

- 4-Hydroxy-3-methoxymethamphetamine
Retention Time: 2.24 min
- Amphetamine
Retention Time: 3.86 min
- Methamphetamine
Retention Time: 4.71 min
- 3,4-Methylenedioxyamphetamine (MDA)-d5
Retention Time: 4.79 min
- 3,4-Methylenedioxyamphetamine (MDA)
Retention Time: 4.85 min
- 3,4-Methylenedioxyamphetamine-d5 (MDMA)-d5
Retention Time: 5.72 min
- 3,4-Methylenedioxyamphetamine (MDMA)
Retention Time: 5.79 min

