

Biogenic Amines in Fish and Fish Products According to § 35 LMBG**Column:** Gemini® 5 µm NX-C18 110 Å, LC Column 250 x 4.6 mm, Ea**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4454-E0**Elution Type:** Gradient**Eluent A:** Water + 100 mM NaOAc, pH 4.5 + 10 mM Sodium 1-octanesulfonate**Eluent B:** Water + 155 mM NaOAc, pH 4.5 + 10 mM Sodium 1-octanesulfonate / ACN (77:23)

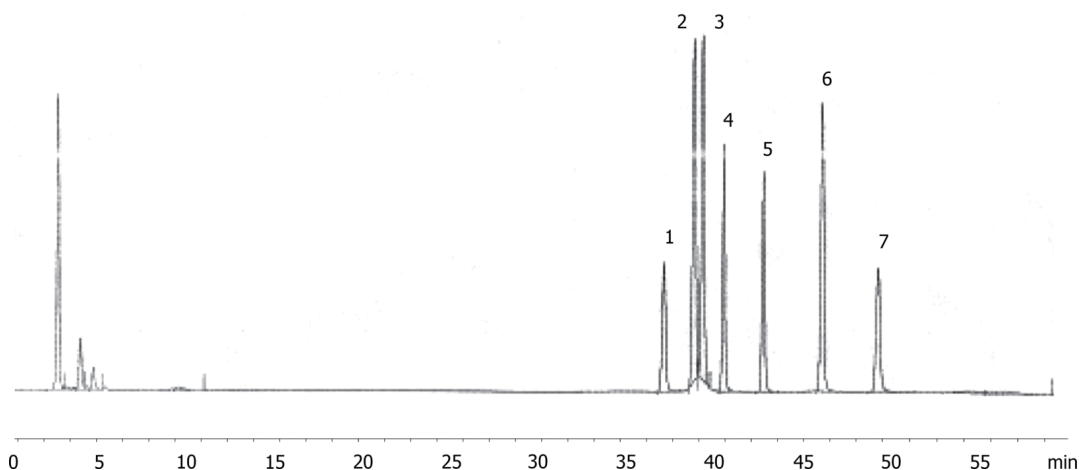
Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	80	20
	2	20	60	40
	3	28	25	75
	4	34	25	75
	5	34.1	0	100
	6	50	0	100
	7	50.1	80	20
	8	60	80	20

Flow Rate: 1 mL/min**Col. Temp.:** 20 °C**Detection:** Fluorescence (FLUOR) @ 330 nm (20 °C)**Detector Info:** Online post-column derivatization with OPA**Analyst Note:** The application was kindly provided by the Amt für Verbraucherschutz in Mettmann, Germany.

17911



Products used in this application:

**ANALYTES:**

- 1 Tyramine
- 2 Putrescine
- 3 Cadaverine
- 4 Histamine
- 5 Agmatine
- 6 Spermidine
- 7 Spermine

