


Hormones/Steroids by MSD on a ZB-5MSi**Column:** ZebtronTM ZB-5MSi, GC Cap. Column 10 m x 0.18 mm x 0.18 μ m, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 10 meters x 0.18 mm x 0.18 μ m**Order No:** 7CD-G018-08**Oven Profile:** 200 °C to 280 °C @ 10 °C/min**Carrier Gas:** Constant Flow Helium, 0.7 mL/min**Injection:** Split 50:1 0.5 μ L @ 260°C**Detection:** Mass Selective (MSD) (320°C)**Analyst Note:** Analytes are 100ppm oxime-TMS derivatized**Products used in this application:****ANALYTES:**

- 1 Androsterone
- 2 17-alpha-Estradiol
- 3 DHEA
- 4 17-beta-Estradiol
- 5 Testosterone
- 6 Testosterone
- 7 Estrone

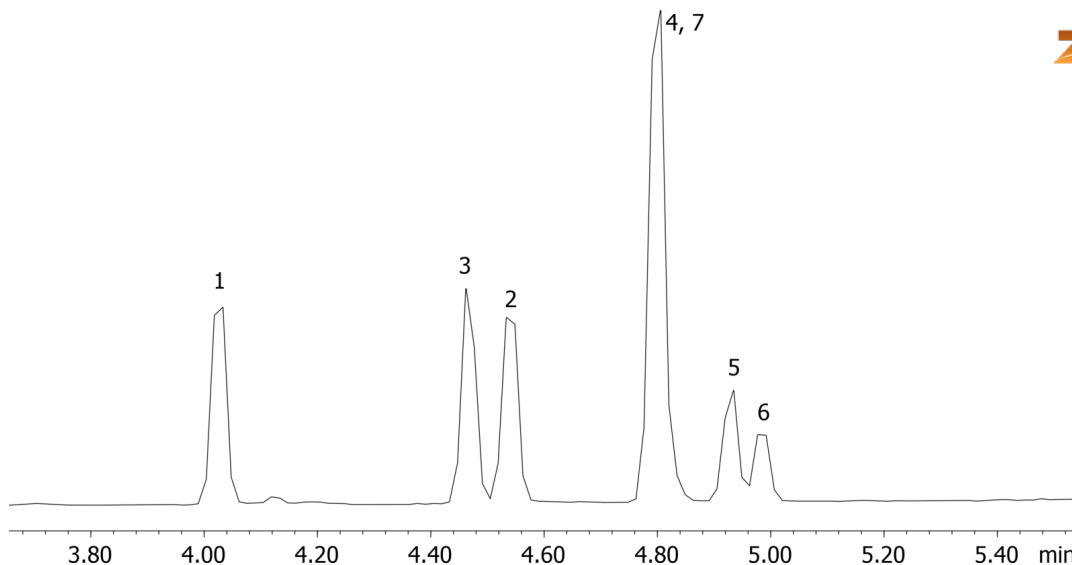