

ABN Mix on Synergi 4u Hydro-RP 150x4.6


Column: Synergi™ 4 µm Hydro-RP 80 Å, LC Column 150 x 4.6 mm, Ea**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4375-E0**Elution Type:** Gradient**Eluent A:** Phosphate Buffer, pH 2.5**Eluent B:** Acetonitrile**Eluent C:** 0**Eluent D:** 0

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B	Pct C	Pct D
	1	0	95	5	0	0
	2	8	20	80	0	0
	3	10	20	80	0	0
	4	10.5	95	5	0	0
	5	12	95	5	0	0

Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Diode Array (PDA) @ 254 210 (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, AQ C18 4 x 3.0mm, 10/Pk Part No.: AJ0-7511

- Holder Part No.: KJ0-4282

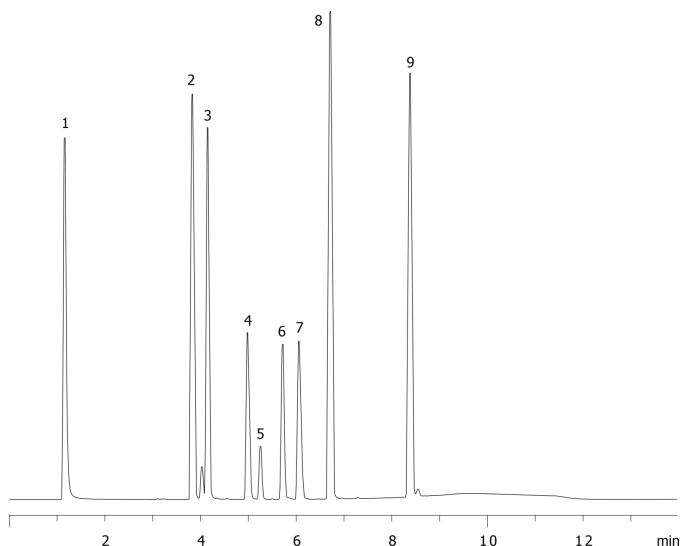

ANALYTES:

- 1 Pyridine
Retention Time: 1.16 min
- 2 Quinidine
Retention Time: 3.83 min
- 3 Sulfathiazole
Retention Time: 4.15 min
- 4 Triprolidine
Retention Time: 4.99 min
- 5 Benzyl alcohol
Retention Time: 5.26 min
- 6 Nortryptiline
Retention Time: 5.72 min
- 7 Phenol
Retention Time: 5.72 min

©2025 Phenomenex Inc. All rights reserved.

For more information contact your Phenomenex Representative at info@phenomenex.com


 3-Methyl 4-benzoic acid
 Phenomenex products are available worldwide.
 Retention Time: 6.71 min
www.phenomenex.com.cn

- 9 3-Methyl salicylaldehyde
Retention Time: 8.32 min

info@phenomenex.com