

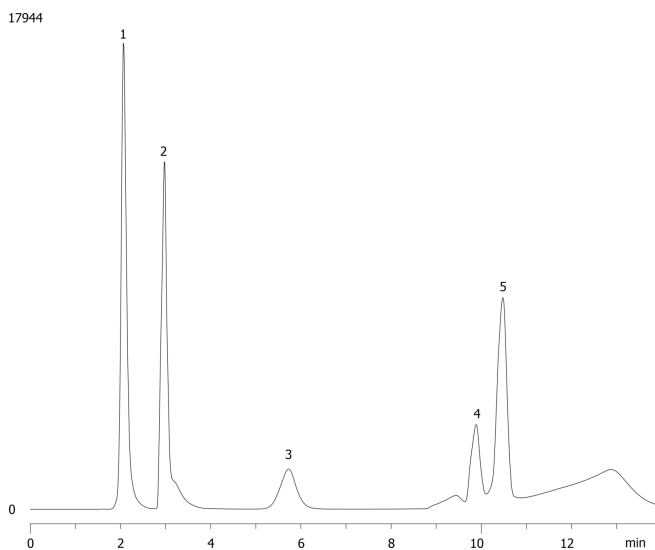
Polar Probes on Synergi Polar**Column:** SynergiTM 4 μ m Polar-RP 80 Å, LC Column 150 x 4.6 mm, Ea**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4336-E0**Elution Type:** Gradient**Eluent A:** Phosphate Buffer, pH 2.5**Eluent B:** Acetonitrile**Eluent C:** 0**Eluent D:** 0

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B	Pct C	Pct D
	1	0	100	0	0	0
	2	5	100	0	0	0
	3	5.01	95	5	0	0
	4	10	70	30	0	0
	5	10.5	100	0	0	0
	6	14	100	0	0	0

Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Diode Array (PDA) @ 210 nm (ambient)**Analyst Note:** SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Polar-RP 4 x 3.0mm ID, 10/Pk Part No.: AJ0-6076

- Holder Part No.: KJ0-4282

**ANALYTES:**

- 1 Cytosine
- 2 Adenine
- 3 Thymine
- 4 Morphine
- 5 Hydromorphone

