

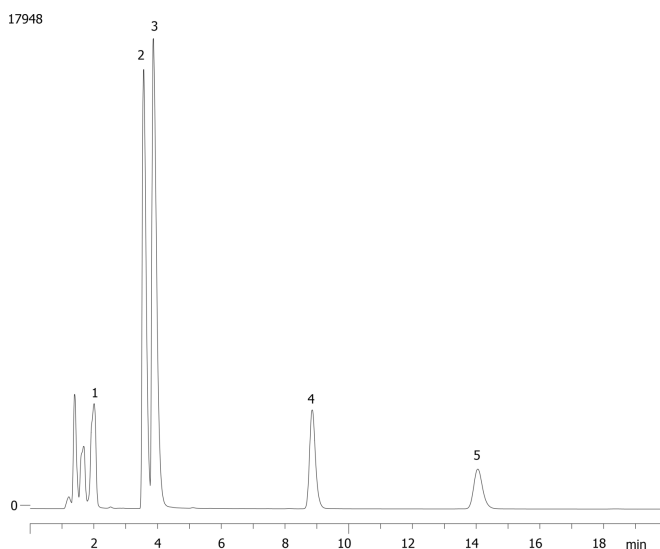
Polar and Non Polar Bases and Hydrophobics on Synergi Fusion**Column:** SynergiTM 4 μ m Fusion-RP 80 Å, LC Column 150 x 4.6 mm, Ea**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4424-E0**Elution Type:** Gradient**Eluent A:** Phosphate Buffer, pH 2.5**Eluent B:** Acetonitrile**Eluent C:** 0**Eluent D:** 0

Gradient	Step No.	Time (min)	Pct A	Pct B	Pct C	Pct D
Profile:	1	20	60	40	0	0

Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Diode Array (PDA) @ 210 nm (ambient)**Analyst Note:** SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Fusion-RP 4 x 3.0mm ID, 10/Pk Part No.: AJ0-7557

- Holder Part No.: KJ0-4282

**ANALYTES:**

- 1 Ranitidine
- 2 Clonidine
- 3 Nortryptiline
- 4 Amytryptiline
- 5 Propyl paraben
- 6 Butyl paraben

