


Blood Alcohols on ZB-BAC-2**Column:** Zebron[™] ZB-BAC-2, GC Cap. Column 30 m x 0.53 mm x 2.00 μ m, Ea**Phase:** Proprietary Blood Alcohol Phase**Dimensions:** 30 meters x 0.53 mm x 2 μ m**Order No:** 7HK-G022-32**Oven Profile:** 40 °C Isothermal**Carrier Gas:** Constant Flow Helium, 80 cm/sec**Injection:** Split 0.8:1 1 μ L @ 150°C**Detection:** Flame Ionization (FID) (250°C)**Analyst Note:** Analytes 0.025% and internal standards 0.100% in water.**Products used in this application:**

17983

**ANALYTES:**

- 1 Methanol
- 2 Acetaldehyde
- 3 Ethanol
- 4 Isopropanol
- 5 Acetone
- 6 tert-Butanol
- 7 n-Propanol
- 8 2-Butanol

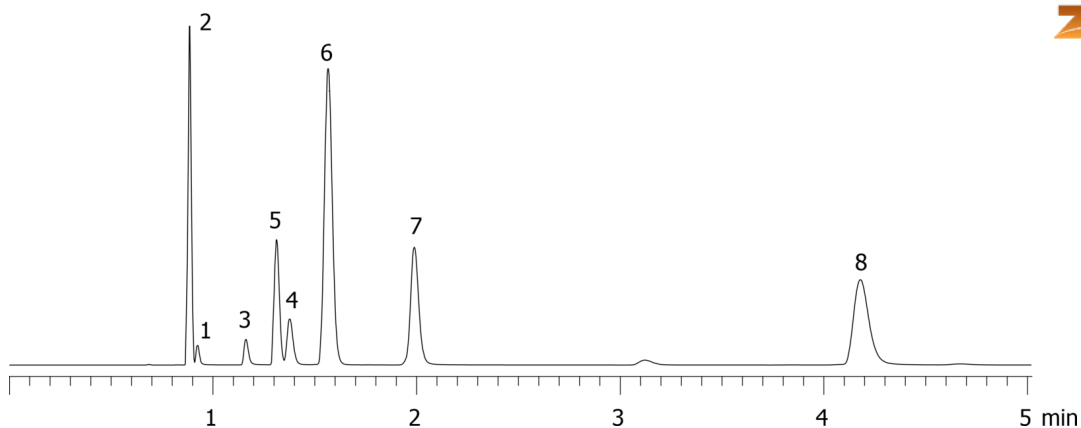