

Nitrophenols in ACN on Synergi 4 μ m Polar-RP 150 x 4.6**Column:** SynergiTM 4 μ m Polar-RP 80 Å, LC Column 150 x 4.6 mm, Ea**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4336-E0**Elution Type:** Isocratic**Eluent A:** Phosphate Buffer**Eluent B:** Acetonitrile

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	72	28

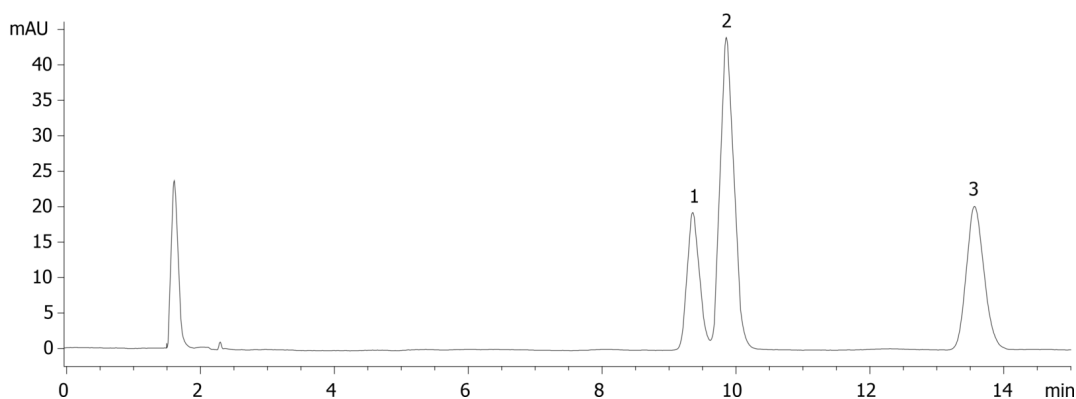
Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ nm (ambient)**Analyst Note:** SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Polar-RP 4 x 3.0mm, 10/Pk Part No.: AJ0-6076

- Holder Part No.: KJ0-4282



17996

**ANALYTES:****1** 4-Nitrophenol

Retention Time: 9.35 min

2 3-Nitrophenol

Retention Time: 9.85 min

3 2-Nitrophenol

Retention Time: 13.56 min

