

Nitrophenols in MeOH on Synergi 4 μ m Polar-RP 150 x 4.6

Column: SynergiTM 4 μ m Polar-RP 80 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4336-E0

Elution Type: Isocratic

Eluent A: Phosphate Buffer

Eluent B: Methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	60	40

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ nm (ambient)


Analyst Note: SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Polar-RP 4 x 3.0mm, 10/Pk Part No.: AJ0-6076

- Holder Part No.: KJ0-4282



17997

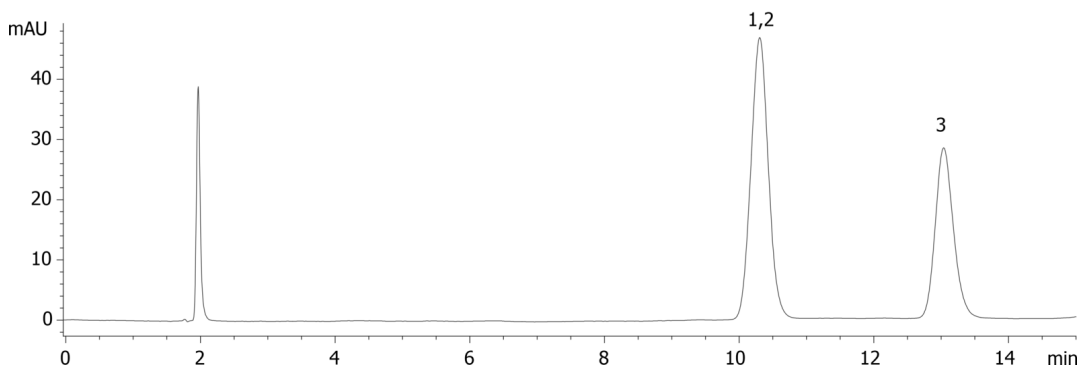

ANALYTES:

1 4-Nitrophenol

Retention Time: 9.99 min

2 3-Nitrophenol

Retention Time: 10.18 min

3 2-Nitrophenol

Retention Time: 13.04 min

