

Prodigy Test Mix in MeOH on Gemini-NX 5 µm 150 x 4.6

Column: Gemini® 5 µm NX-C18 110 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4454-E0

Elution Type: Isocratic

Eluent A: Water

Eluent B: Methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	25	75

Flow Rate: 1.5 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Gemini-NX C18 4 x 3.0mm, 10/Pk Part No.: AJ0-8368


- Holder Part No.: KJ0-4282



Products used in this application:



17999

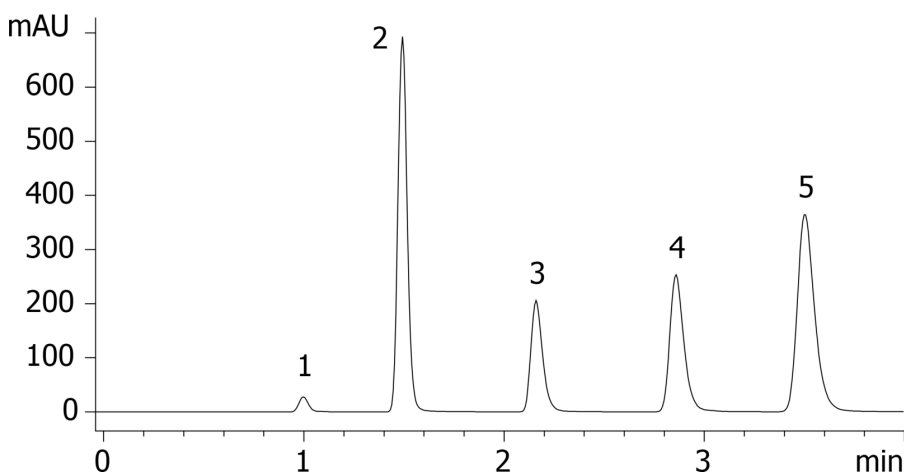

ANALYTES:

- 1** Uracil
Retention Time: 1 min
- 2** Acetophenone
Retention Time: 1.49 min
- 3** Benzene
Retention Time: 2.16 min
- 4** Toluene
Retention Time: 2.86 min
- 5** Naphthalene
Retention Time: 3.5 min

