


Blood Alcohols using Hydrogen on ZB-BAC-1**Column:** Zebron[™] ZB-BAC-1, GC Cap. Column 30 m x 0.53 mm x 3.00 μ m, Ea**Phase:** Proprietary Blood Alcohol Phase**Dimensions:** 30 meters x 0.53 mm x 3 μ m**Order No:** 7HK-G021-36**Oven Profile:** 40 °C Isothermal**Carrier Gas:** Constant Flow Hydrogen, 80 cm/sec**Injection:** Split 5:1 1 μ L @ 150°C**Detection:** Flame Ionization (FID) (250°C)**Analyst Note:** All analytes are 0.100% in water**Products used in this application:**

18025

**ANALYTES:**

- 1 Methanol
- 2 Acetaldehyde
- 3 Ethanol
- 4 Isopropanol
- 5 Acetone
- 6 n-Propanol

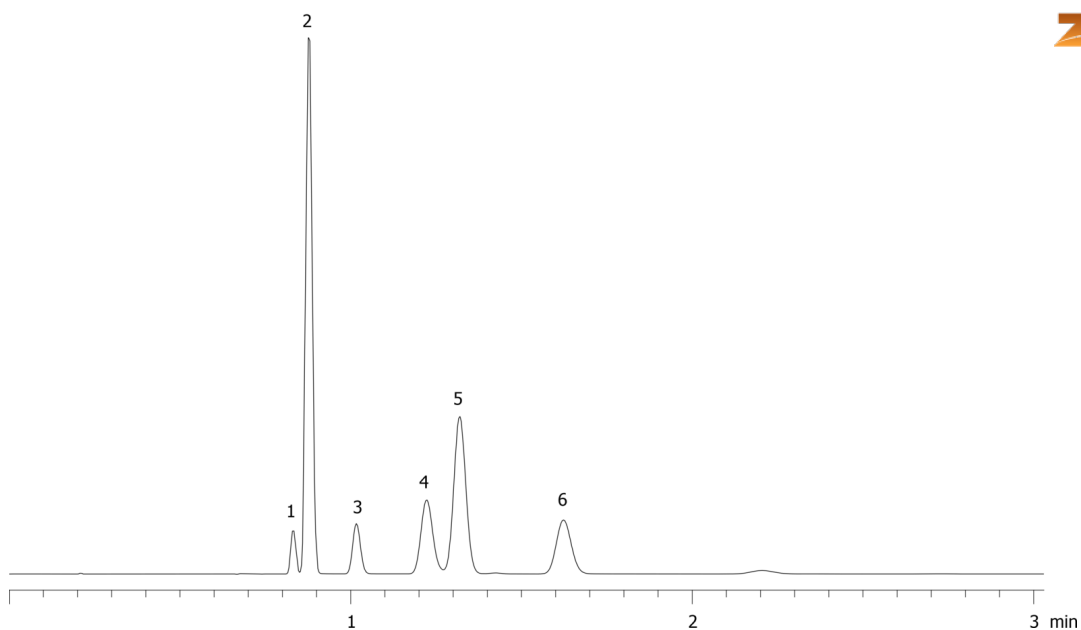