

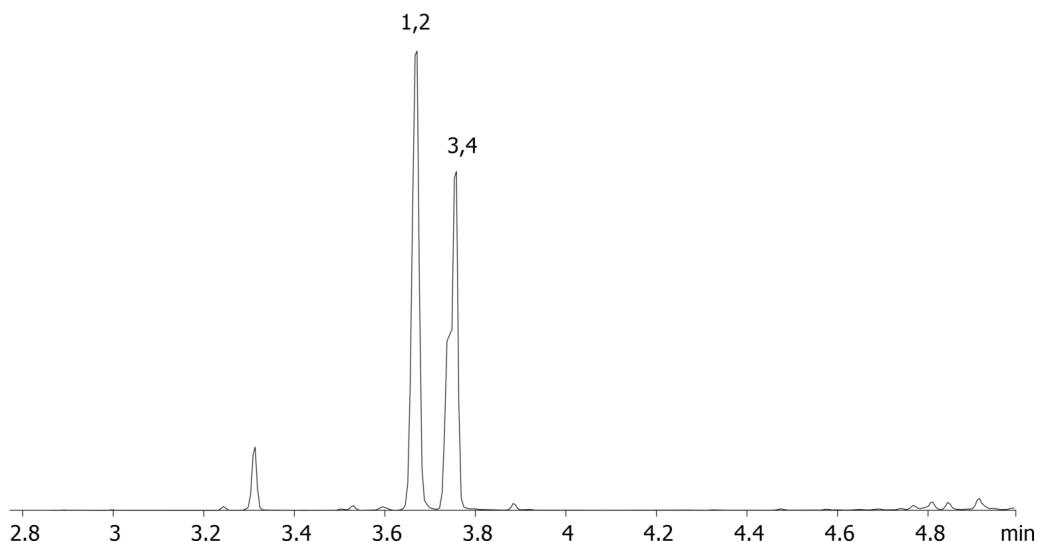
Opiates from urine by Screen C on a ZB-Drug1**Column:** Zebtron[™] ZB-Drug-1, GC Cap. Column 10 m x 0.18 mm x 0.18 µm, Ea**Phase:** Proprietary Drug Phase**Dimensions:** 10 meters x 0.18 mm x 0.18 µm**Order No:** 7CD-G023-08**Oven Profile:** 180 °C to 320 °C @ 25 °C/min for 1.4 min**Carrier Gas:** Constant Flow Helium, 0.7 mL/min**Injection:** Split 10:1 1 µL @ 280°C**Detection:** Mass Selective (MSD) (320°C)**Analyst Note:** Analytes are 2000 ng/mL in urine derivatized with BSTFA

Sim ions:

(Quant ion*) D6-Codeine (*377,349): Codeine (*371,343,234): D3-Morphine (*432,417): Morphine (*429,414,401)

**Products used in this application:**

18035

**ANALYTES:**

- 1 Codeine-d6
- 2 Codeine
- 3 Morphine-d3
- 4 Morphine

