

Kinetex 2.6µm C18 50 x 2.1mm ID mm van Deemter Curve


Column: Kinetex® 2.6 µm C18 100 Å, LC Column 50 x 2.1 mm, Ea
Dimensions: 50 x 2.1 mm ID
Order No: 00B-4462-AN
Elution Type: Isocratic
Eluent A: Water
Eluent B: Acetonitrile
Gradient Profile:

Step No.	Time (min)	Pct A	Pct B
1	0	50	50

Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ nm (ambient)
18117



Products used in this application:

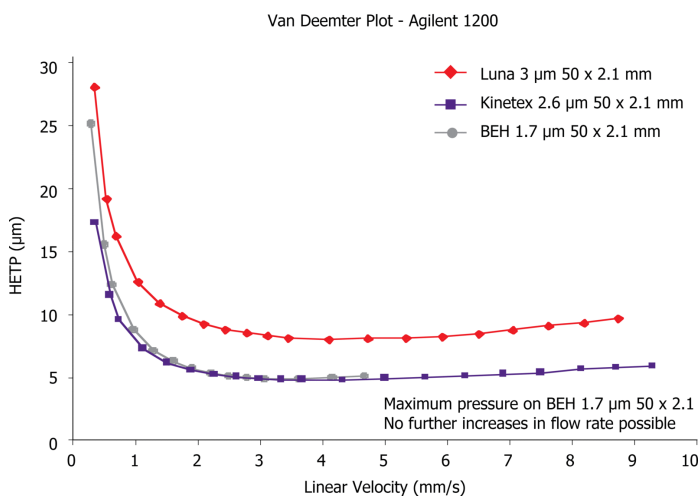

ANALYTES:

- 1 Uracil
- 2 Acetophenone
- 3 Benzene
- 4 Toluene
- 5 Naphthalene

