

Oxazepam on Lux 5 µm Cellulose-1 250 x 4.6mm ID in PO-95%CH₃CN**Column:** Lux® 5 µm Cellulose-1, LC Column 250 x 4.6 mm, Ea**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4459-E0**Elution Type:** Isocratic**Eluent A:** 0.1 % DEA in CH₃CN**Eluent B:** 0.1% DEA in IPA

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	10	95	5

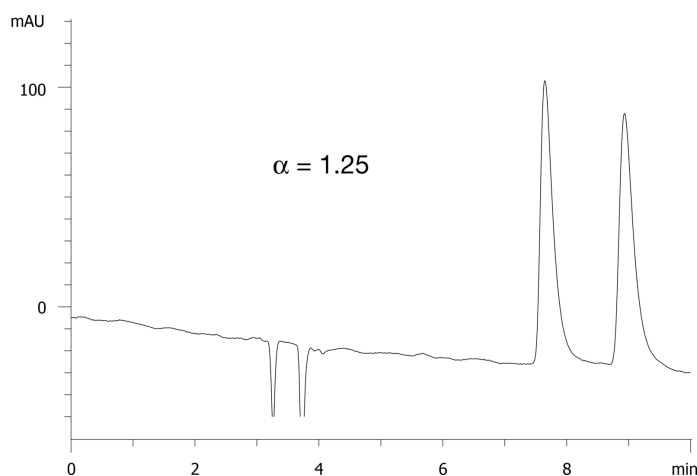
Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Cellulose-1, 4 x 3.0mm ID, 10/Pk Part No.: AJ0-8403

- Holder Part No.: KJ0-4282

**Products used in this application:**

18141

**ANALYTES:****1** Oxazepam

Retention Time: 7.65 min

2 Oxazepam

Retention Time: 8.93 min

