


Aroclor 1242: DIN 51527 on ZB-35HT**Column:** Zebron[™] ZB-35HT, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** 35% Phenyl 65% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G025-11**Oven Profile:** 50 °C for 0.5 min to 220 °C @ 30 °C/min to 320 °C @ 6 °C/min for 2 min**Carrier Gas:** Constant Flow Helium, 0.5 mL/min**Injection:** Splitless :1 1 μ L @ 250°C**Detection:** Electron Capture (ECD) (350°C)**Detector Info:** Used a pressure pulse at 40 psi for 0.25 min**Analyst Note:** Total concentration of Aroclors was 90 ppm in iso-octane**Products used in this application:**

Aroclor 1242: DIN 51527 on ZB-35HT

ANALYTES:

- 1 TCMX
- 2 BZ# 31
- 3 BZ# 28
- 4 BZ# 52
- 5 BZ# 101
- 6 BZ# 77
- 7 BZ# 123
- 8 BZ# 118
- 9 BZ# 153
- 10 BZ# 114
- 11 BZ# 105
- 12 BZ# 138
- 13 BZ# 126
- 14 BZ# 167
- 15 BZ# 156
- 16 BZ# 180
- 17 BZ# 157
- 18 BZ# 169
- 19 BZ# 170
- 20 BZ# 189
- 21 DCB

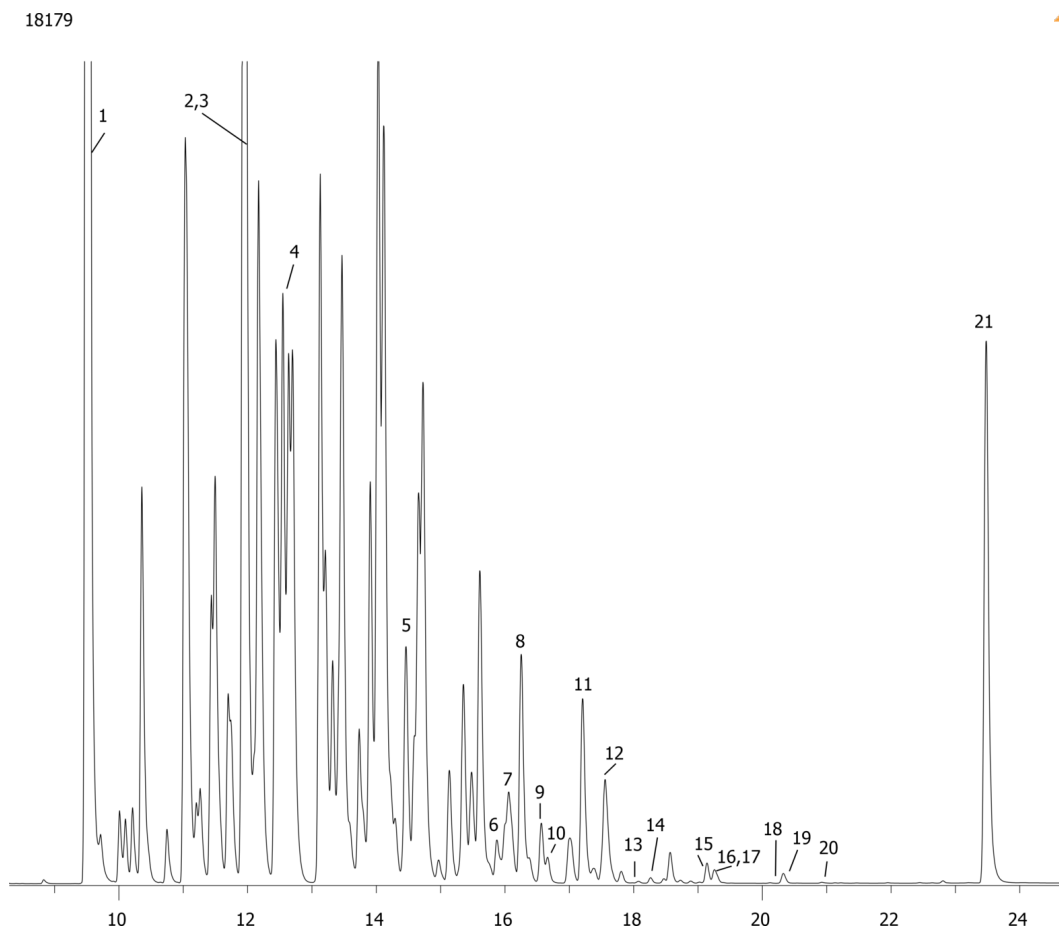