

Bendroflumethiazide on Lux 5u Cellulose-2 by SFC**Column:** Lux® 5 µm Cellulose-2, LC Column 250 x 4.6 mm, Ea**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4457-E0**Elution Type:** Isocratic**Eluent A:** 0.1% DEA+0.1% TFA in methanol**Eluent B:** CO₂

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	30	70

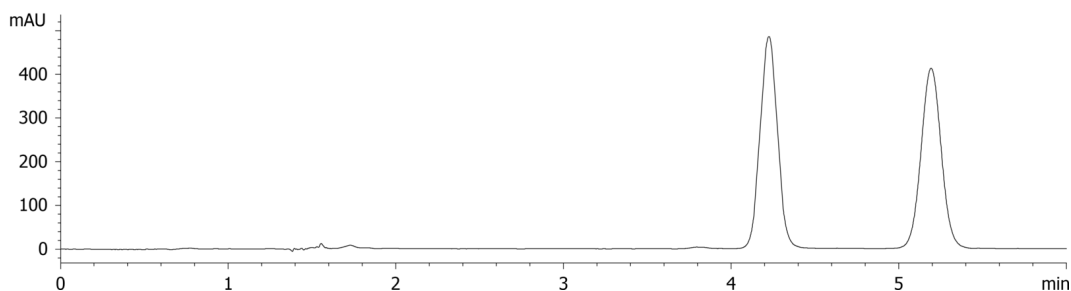
Flow Rate: 2 mL/min**Col. Temp.:** 40 °C**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 273 nm (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Cellulose-2, 4 x 3.0mm ID, 10/Pk Part No.: AJ0-8366

- Holder Part No.: KJ0-4282

**Products used in this application:**

18207

**ANALYTES:****1** Bendroflumethiazide