

Trans Stilbene Oxide on a Lux 5u Cellulose 2 Axia 50mm ID column

Column: Lux® 5 µm Cellulose-2, LC Column 250 x 50 mm, AXIA™ Packed, Ea

Dimensions: 250 x 50 mm ID

Order No: 00G-4457-V0-AX

Elution Type: Isocratic

Eluent A: Hexane

Eluent B: IPA

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	1	0	90	10


Flow Rate: 50 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)
18212



Products used in this application:

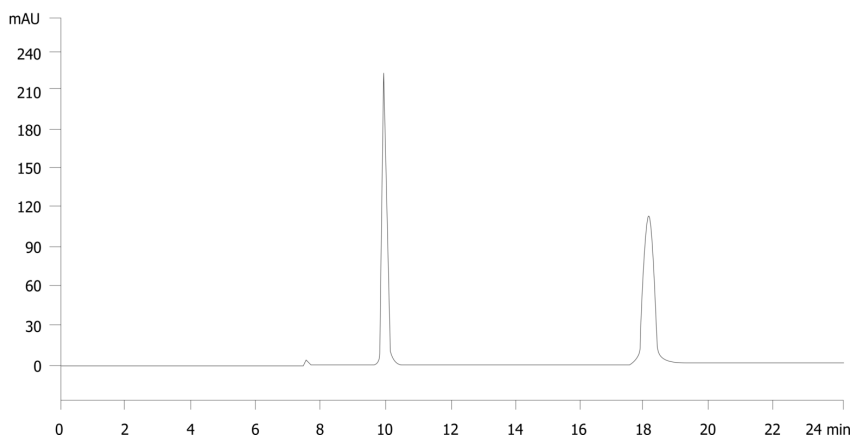

ANALYTES:

1 trans-Stilbene oxide

