

Trans Stilbene Oxide on a Lux 5u Cellulose 2 Axia 30mm ID column

Column: Lux® 5 µm Cellulose-2, LC Column 250 x 30 mm, AXIA™ Packed, Ea

Dimensions: 250 x 30 mm ID

Order No: 00G-4457-U0-AX

Elution Type: Isocratic

Eluent A: Hexane

Eluent B: IPA

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	0	90	10	

Flow Rate: 40 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard PREP Cartridge, Lux Cellulose-2, 15 x 30mm ID, Ea Part No.: AJ0-8401


- Holder Part No.: AJ0-8277



Products used in this application:



18213

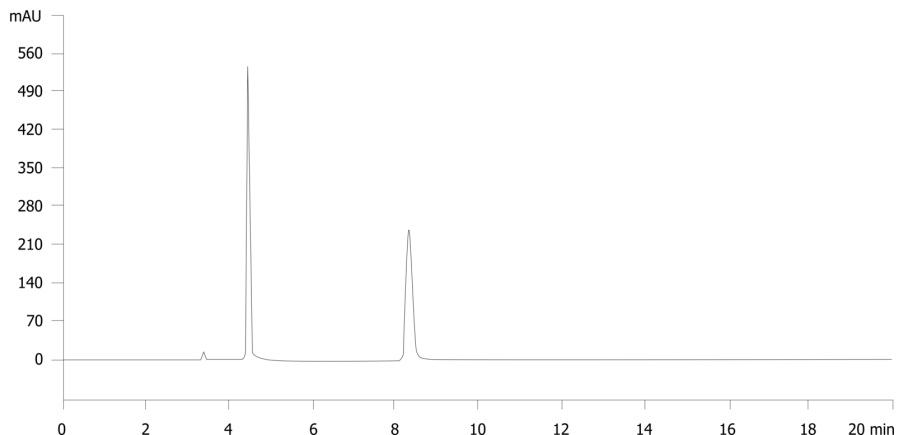

ANALYTES:

1 trans-Stilbene oxide

