

Troger's Base on Lux 3u Cellulose-1

Column: Lux® 3 µm Cellulose-1, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4458-E0

Elution Type: Isocratic

Eluent A: Hexane

Eluent B: Isopropanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	0	90	10	

Flow Rate: 1 mL/min

Col. Temp.: 25 °C

Detection: UV-Vis Abs.-Diode Array (PDA) @ 230 nm (ambient)

Detector Info: alpha=1.24 and Rs=5.22

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Cellulose-1, 4 x 3.0mm, 10/Pk Part No.: AJ0-8403


- Holder Part No.: KJ0-4282



Products used in this application:



18216

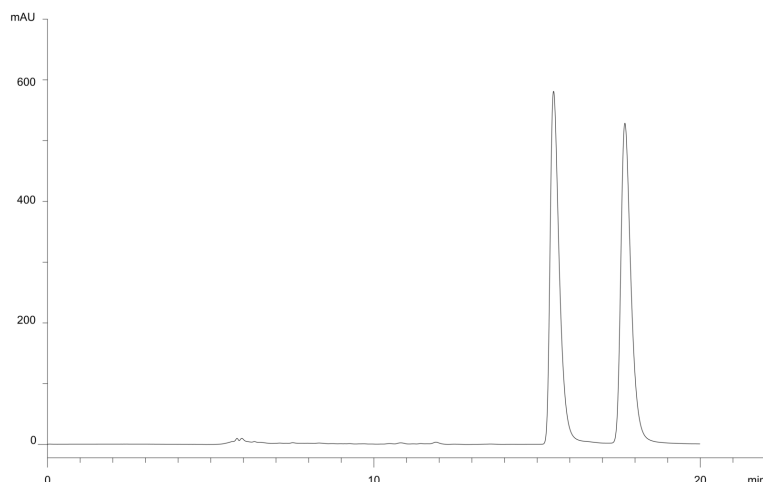

ANALYTES:

1 Troger's base

