

HPLC Application

ID No.: 18220

Sensitivity - Kinetex 2.6 μ m C18 50 x 2.1mm ID (0.5 mL/min)

Column: Kinetex[®] 2.6 μ m C18 100 Å, LC Column 50 x 2.1 mm, Ea

Dimensions: 50 x 2.1 mm ID

Order No: 00B-4462-AN

Elution Type: Isocratic

Eluent A: Water

Eluent B: Acetonitrile


Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	50	50

Flow Rate: 0.5 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (25 °C)

18220

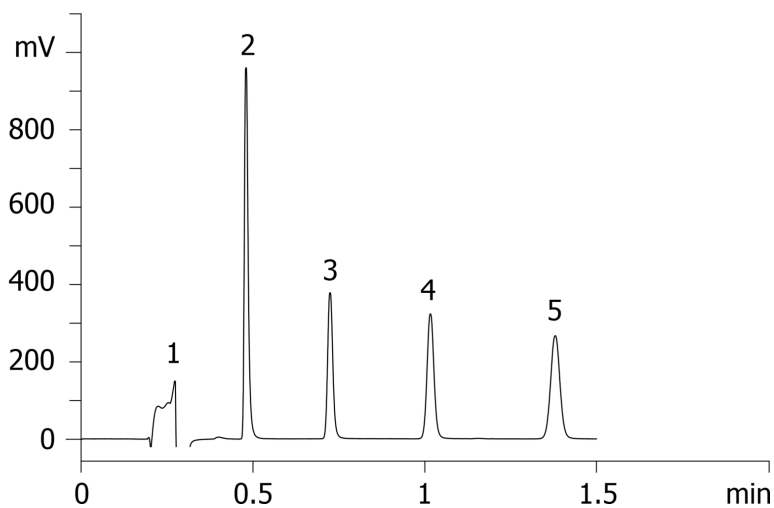

ANALYTES:

- 1 Uracil
- 2 Acetophenone
- 3 Benzene
- 4 Toluene
- 5 Naphthalene



Products used in this application:

