

Pharma ABN Test 12 Mix on Luna 5um C18(2) 250x4.6mm 1mL/min

Column: Luna® 5 µm C18(2) 100 Å, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4252-E0

Elution Type: Gradient

Eluent A: water w/ 0.1% formic acid

Eluent B: acetonitrile w/ 0.1% formic acid

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	2.79	95	5
	3	36.15	5	95
	4	36.38	95	5
	5	50.93	95	5

Flow Rate: 1 mL/min

Col. Temp.: 45 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (25 °C)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287


- Holder Part No.: KJ0-4282



Products used in this application:



18230

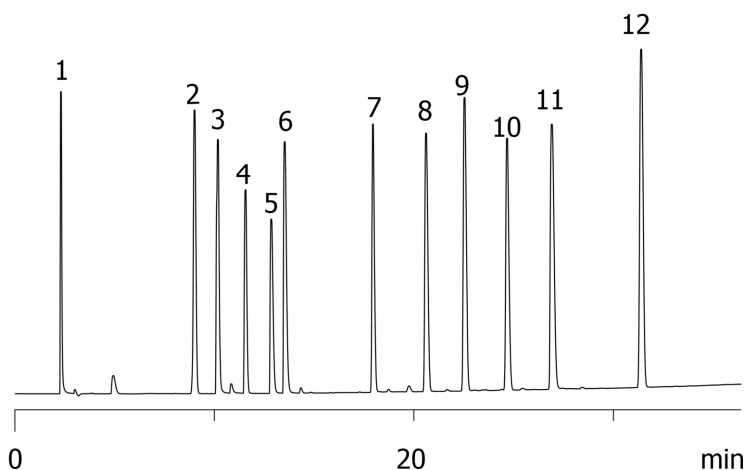

Pharma ABN Test 12 Mix on Luna 5um C18(2) 250x4.6mm 1mL/min

ANALYTES:

- 1** Pyridine
- 2** Acetaminophen
- 3** Quinine
- 4** Acebutolol
- 5** Chlorpheniramine
- 6** Triprolidine
- 7** Prednisolone
- 8** 4-Chlorobenzoic acid
- 9** 4-Chlorocinnamic acid
- 10** Diazepam
- 11** Diflunisal
- 12** Hexanophenone

