

Sulfa Drugs on Kinetex 1.7µm C18 50 x 2.1mm ID (0.6 mL/min)

Column: Kinetex® 1.7 µm C18 100 Å, LC Column 50 x 2.1 mm, Ea

Dimensions: 50 x 2.1 mm ID

Order No: 00B-4475-AN

Elution Type: Gradient

Eluent A: water with 0.1% formic acid

Eluent B: acetonitrile with 0.1% formic acid

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	2.33	5	95
	3	2.36	95	5
	4	3.37	95	5

Flow Rate: 0.6 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 260 nm (25 °C)

18233



Products used in this application:

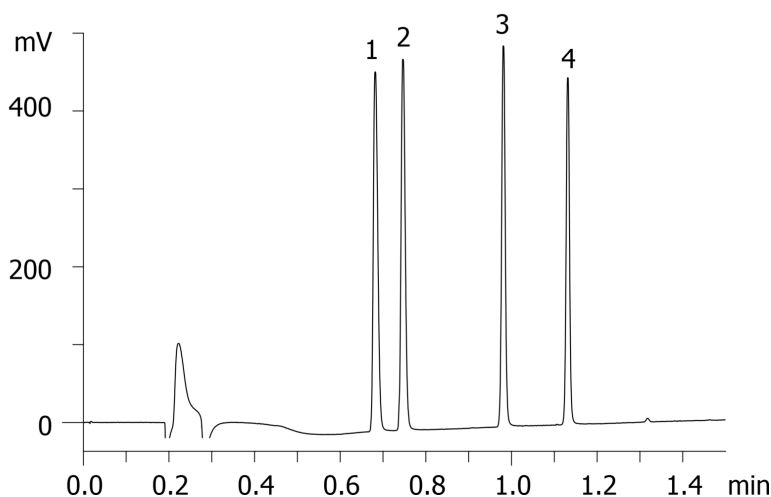

ANALYTES:

- 1 Sulfathiazole
- 2 Sulfamerazine
- 3 Sulfamethoxazole
- 4 Sulfaquinoxaline

