

Impurity Profile on Kinetex 2.6µm C18 150 x 4.6mm ID (1.5mL/min)

Column: Kinetex® 2.6 µm C18 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4462-E0

Elution Type: Gradient

Eluent A: water

Eluent B: acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	1.16	95	5
	3	14.51	5	95
	4	14.6	5	95
	5	14.61	95	5
	6	20.42	95	5

Flow Rate: 1.5 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (25 °C)

18262



Products used in this application:

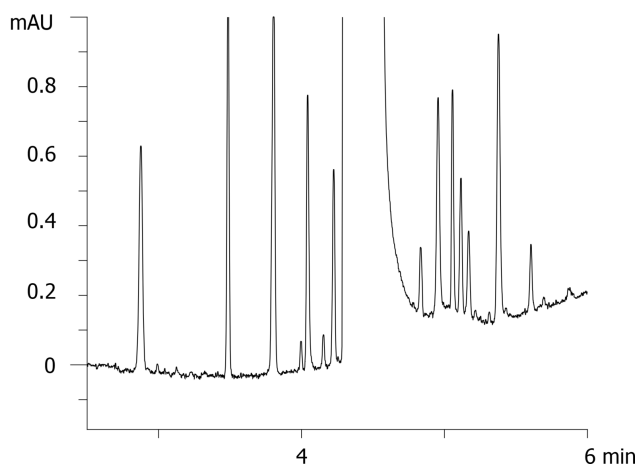

Impurity Profile on Kinetex 2.6µm C18 150 x 4.6mm ID (1.5mL/min)

ANALYTES:

- 1** Pharmaceutical impurities
- 2** Acetaminophen
- 3** Pindolol
- 4** Quinine
- 5** Acebutolol
- 6** Chlorpheniramine
- 7** Triprolidine
- 8** Prednisolone
- 9** Nortriptyline
- 10** 4-Chlorobenzoic acid
- 11** 5-Methyl-2-hydroxy benzaldehyde
- 12** 4-Chlorocinnamic acid
- 13** Diazepam
- 14** Diflunisal
- 15** Niflumic acid
- 16** Hexanophenone

