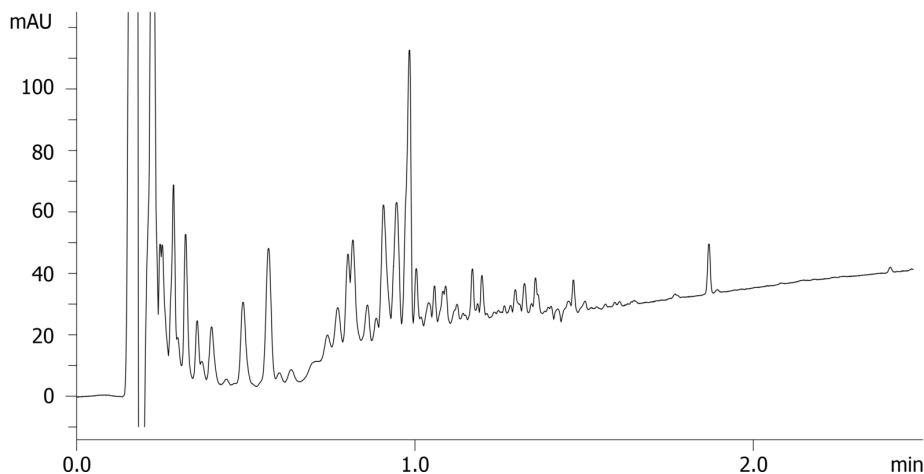


Metabolomics on Kinetex 2.6µm C18 50 x 4.6mm ID**Column:** Kinetex® 2.6 µm C18 100 Å, LC Column 50 x 4.6 mm, Ea**Dimensions:** 50 x 4.6 mm ID**Order No:** 00B-4462-E0**Elution Type:** Gradient**Eluent A:** water with 0.1% formic acid**Eluent B:** acetonitrile with 0.1% formic acid

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	0.2	95	5
	3	2.47	5	95
	4	2.48	95	5
	5	3.47	95	5

Flow Rate: 3 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Diode Array (PDA) @ 220 nm (25 °C)

18275

**Products used in this application:****ANALYTES:****1** Human urine