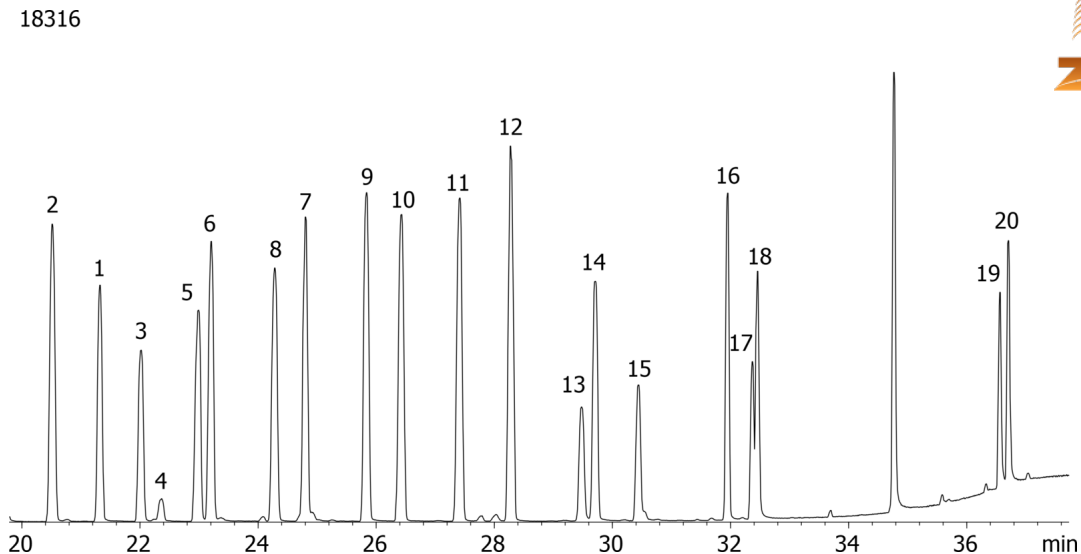


Halogenated Pesticides Mix # 6 on ZB-5MSi

Column: Zebtron[™] ZB-5MSi, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea
Phase: 5% Phenyl 95% Dimethylpolysiloxane
Dimensions: 30 meters x 0.25 mm x 0.25 μ m
Order No: 7HG-G018-11
Oven Profile: 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 5 min
Carrier Gas: Constant Flow Helium, 0.8 mL/min
Injection: On-Column :1 0.5 μ L @ 43°C
Detection: Mass Selective (MSD) (335°C)
Analyst Note: Analytes 10 ppm in Acetone

**Products used in this application:**

Halogenated Pesticides Mix # 6 on ZB-5MSi

ANALYTES:

- 1 Ethalfluralin
- 2 2,3,5,6-Tetrachloroaniline
- 3 Di-allate
- 4 Di-allate
- 5 Simazine
- 6 Atrazine
- 7 Tefluthrin
- 8 Chlordene
- 9 Dimethachlor
- 10 Heptachlor
- 11 Pentachlorothioanisole
- 12 Dacthal
- 13 Pyrifenox
- 14 Fipronil
- 15 Pyrifenox
- 16 Nitrofen
- 17 Etaconazole
- 18 Etaconazole
- 19 Flucythrinate
- 20 Flucythrinate

