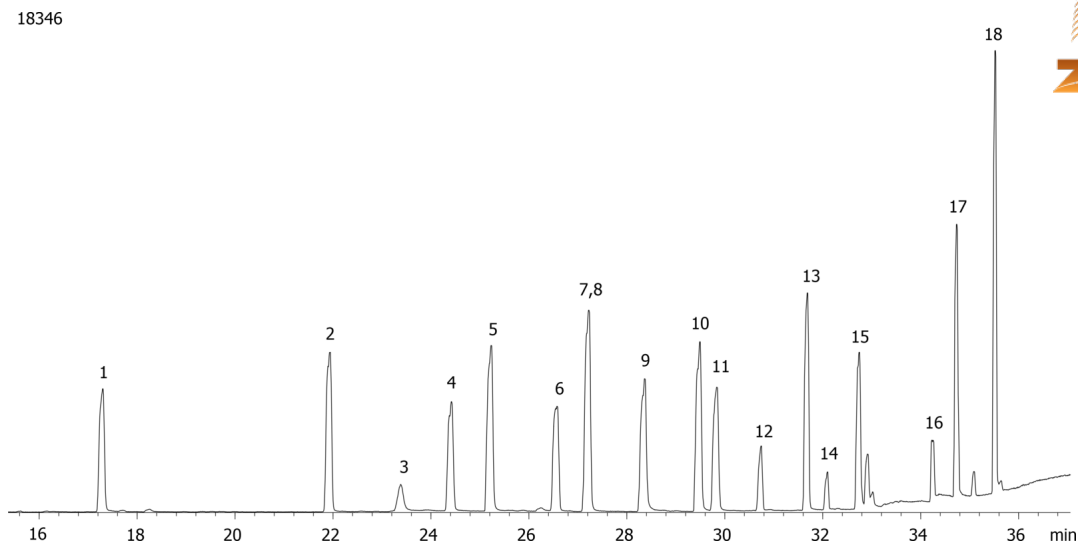


**NSO (Nitrogen/Sulfur/Oxygen) Pesticides Mix # 5 on ZB-XLB**

**Column:** Zebron<sup>™</sup> ZB-XLB, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea  
**Phase:** Proprietary XLB Phase  
**Dimensions:** 30 meters x 0.25 mm x 0.25 µm  
**Order No:** 7HG-G019-11  
**Oven Profile:** 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 5 min  
**Carrier Gas:** Constant Flow Helium, 0.8 mL/min  
**Injection:** On-Column :1 0.5 µL @ 43°C  
**Detection:** Mass Selective (MSD) (335°C)  
**Analyst Note:** Analytes 10 ppm in Acetone

**Products used in this application:**

## NSO (Nitrogen/Sulfur/Oxygen) Pesticides Mix # 5 on ZB-XLB

### ANALYTES:

- 1 Metolcarb
- 2 2,3,5-Trimethacarb
- 3 Desmedipham
- 4 Cycluron
- 5 Secbumetone
- 6 Dioxacarb
- 7 Prometryne
- 8 3-Hydroxycarbofuran
- 9 Naphaleneacetamide
- 10 Tetraconazole
- 11 Pyracarbolid
- 12 Dinobuton
- 13 Napropamide
- 14 Pyrethrin 1
- 15 Tricyclazole
- 16 Nitralin
- 17 Fenoxycarb
- 18 Amitraz

