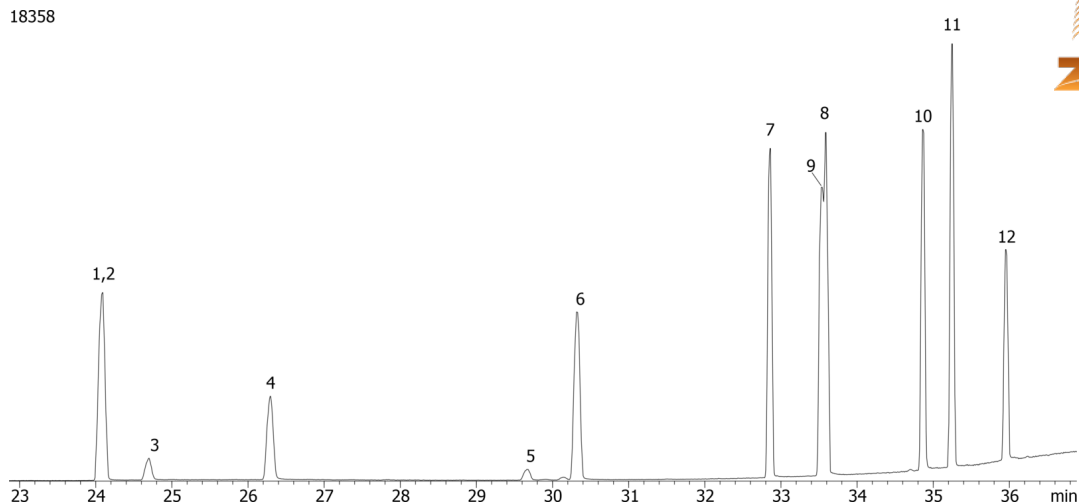


Organophosphorous Pesticides (EPA 8141) Mix # 3 on ZB-XLB

Column: Zebron[™] ZB-XLB, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea
Phase: Proprietary XLB Phase
Dimensions: 30 meters x 0.25 mm x 0.25 µm
Order No: 7HG-G019-11
Oven Profile: 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 5 min
Carrier Gas: Constant Flow Helium, 0.8 mL/min
Injection: On-Column :1 0.5 µL @ 43°C
Detection: Mass Selective (MSD) (335°C)
Analyst Note: Analytes 10 ppm in Acetone

**Products used in this application:**

Organophosphorous Pesticides (EPA 8141) Mix # 3 on ZB-XLB

ANALYTES:

- 1 Terbufos
- 2 Dioxathion
- 3 Phosphamidon
- 4 Phosphamidon
- 5 Chlorfenvinphos
- 6 Chlorfenvinphos
- 7 Ethion
- 8 Carbophenothion
- 9 Famphur
- 10 Phosmet
- 11 Leptophos
- 12 Azinphos-ethyl

